PHILIPS

HF330 I


PHILIPS

English**Page 4**


- Keep page 3 open when reading these operating instructions.

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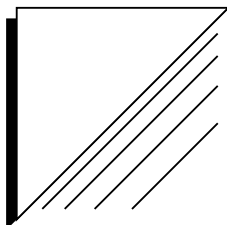
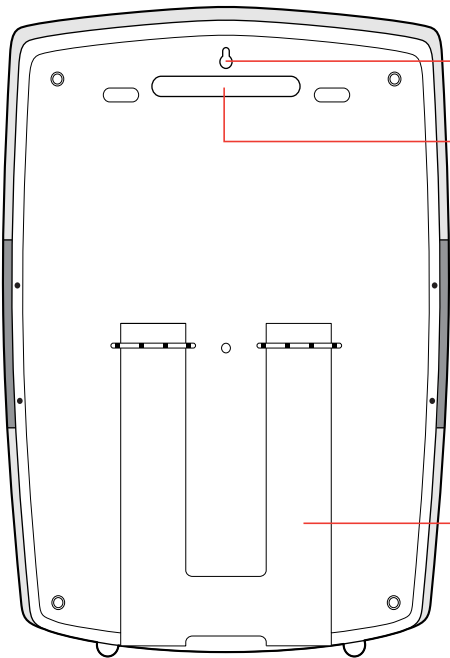
H

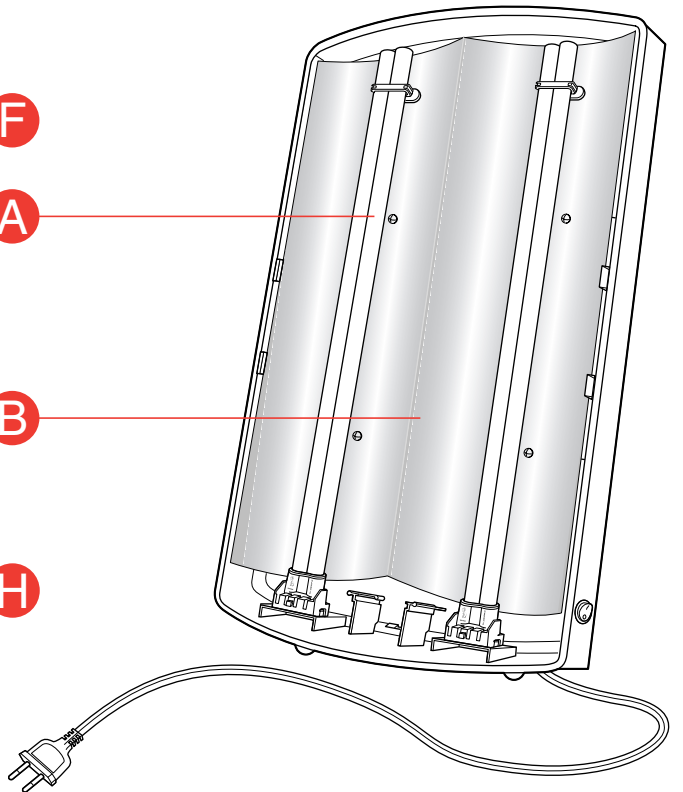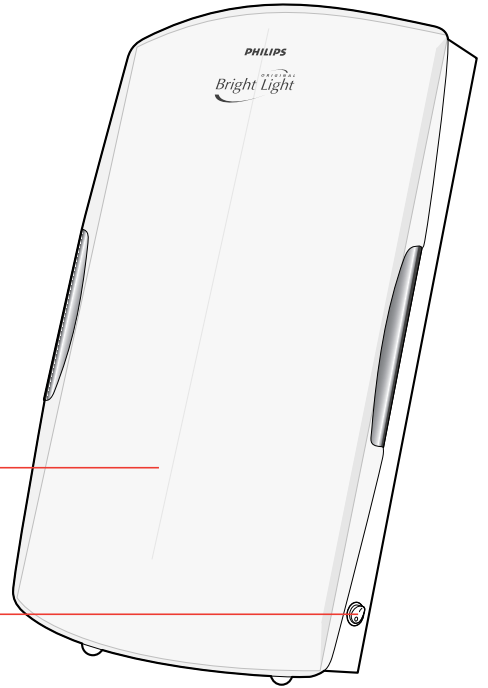Important


Read these instructions for use carefully before using the appliance and save them for future reference.

- Check if the voltage indicated on the appliance corresponds to the local mains voltage before you connect the appliance.
- Keep this appliance away from water! Do not use this appliance in damp environments (e.g. in a bathroom or near a swimming pool).
- If the mains cord is damaged, it must be replaced by Philips, a service centre authorised by Philips or similarly qualified persons in order to avoid a hazard.
- Always unplug the appliance after use.
- Keep the appliance out of the reach of children.
- The duration of the light bath depends on the distance at which the appliance is used. See chapter 'Light bath' for the appropriate distance and duration. Repeat the light bath on at least five successive days (take them at a time that suits you and make this a daily habit during the course). You can start a new light bath course (5 or more successive days) whenever you feel you need it.
- Use the light of this appliance as an addition to the normal light (sunlight, electric light) in your home. Do not use the appliance in a room that is not illuminated by any other light source.
- It is not advisable to look into the light too long.
- After the first few light baths you may experience some eye strain and/or a light headache. These side effects tend to disappear quickly.
- Consult your doctor before you start using the appliance if:
 - you have suffered from or are suffering from severe depression;
 - you are hypersensitive to light (e.g. diabetics and epileptics);
 - you have eye complaints;
 - you are on specific medication (i.e. certain antidepressants, psychotropic drugs or malaria tablets);
 - you have high blood pressure.
- If using the appliance causes you to become too active, we advise you to reduce the light dose. If there is no improvement, we advise you to consult your doctor.
- If you feel exhausted, distressed or restless and/or suffer from sleeping problems (e.g. insomnia) after having used the appliance, it is always advisable to consult your doctor or psychotherapist.
- This appliance meets the requirements of MDD 93/42/EEC for medical appliances.


General description (fig. 1)

- A Tube lamps (type Philips Original Bright Light PL-L 55 W)
- B Reflector
- C Opalescent/satin screen (acrylic)
- D On/Off switch
- E Suspension hole
- F Stand (fold-out)
- G Handle (for transport)
- H Screw and plug (for hanging the appliance)



Making the appliance ready for use

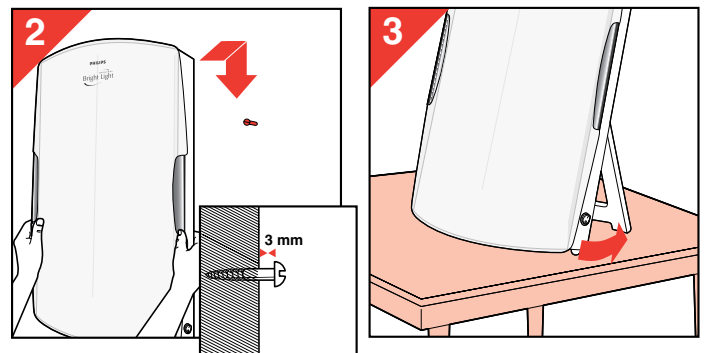
- Hang the appliance on the wall (fig. 2) or unfold the stand (fig. 3) and place it on a table top.
- Insert the plug into a wall socket (220-240V) (fig. 4).

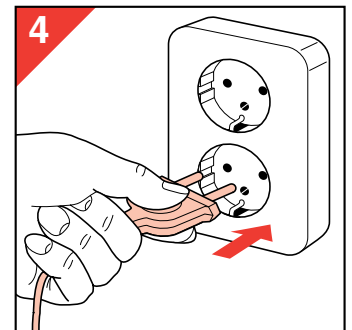Light bath

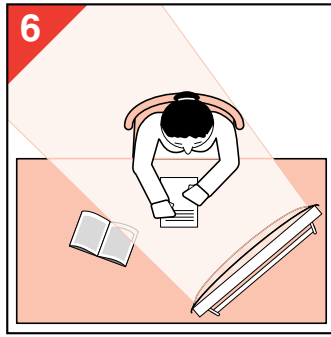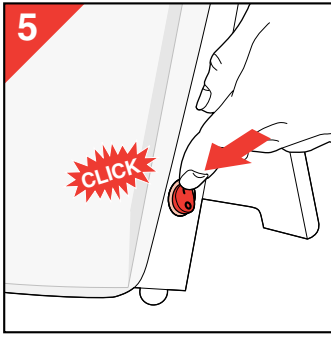- Switch on (fig. 5).
- Position yourself in front of the appliance (fig. 6).
For a light bath, the distance between your face and the opalescent/satin screen should be between 25 and 75 cm.

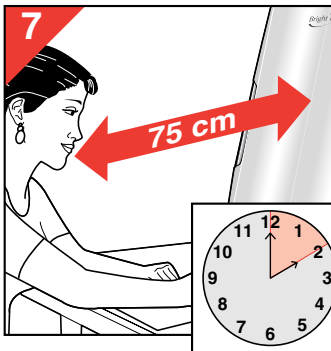- The advisable duration of the lightbath depends on the selected operational distance:
approx. **2 hours** for a distance of approx. **75 cm** (fig. 7);
approx. **1 hour** for a distance of approx. **50 cm** (fig. 8);
approx. **½ hour** for a distance of approx. **25 cm** (fig. 9).


- You need not take the whole lightbath without a break. If it suits your schedule better, you may interrupt your lightbath and continue it later. Even the 30 minute session may be split into two sessions if you wish.
- You may look into the lamps every now and then, not continuously. So you are free to engage in such activities as reading, writing, fancy-work or handiwork, and look into the lamps every now and then (figs. 10, 11).

- Preferably have your lightbaths between 8 o'clock in the morning (a.m.) and 8 o'clock in the evening (p.m.) (fig. 12).



- Repeat the lightbaths until you feel better. Generally, improvements occur within one week. So as a rule, we advise you to repeat the lightbath on at least 5 successive days (fig. 13).

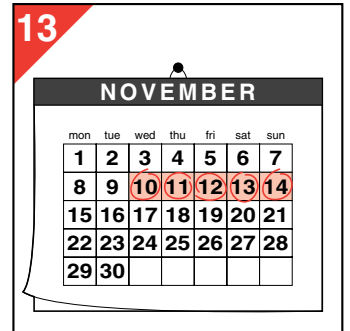- You may repeat the light bath as much as you like, during the 'dark period' of the year. Do repeat the lightbath on at least five successive days.

After use



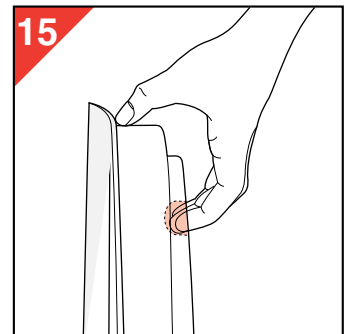
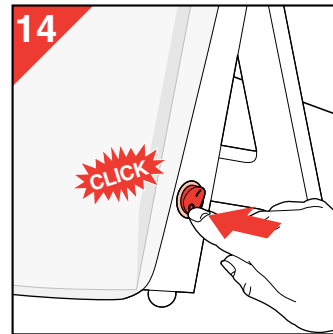
- Switch the appliance off (fig. 14) and remove the plug from the wall socket.
- You can lift the appliance by its handle (fig. 15).

Cleaning

Always unplug the appliance and let it cool down before cleaning.

- You can clean the outside of the appliance with a damp cloth (fig. 16). Make sure that no water runs into the appliance (fig. 17). Do not use abrasives or scourers, or liquids such as alcohol, petrol or acetone.

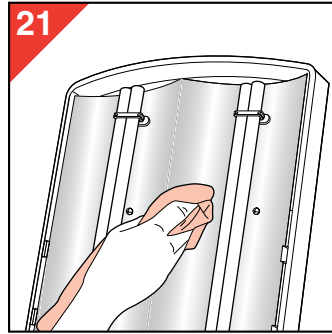
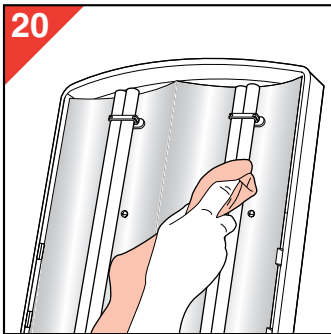
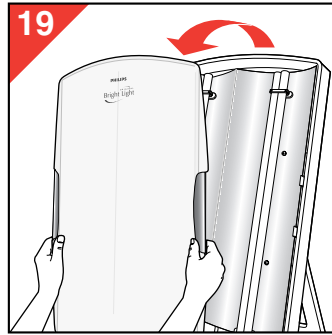
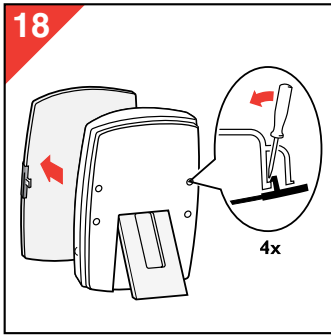
- From time to time the tube lamps, reflectors and the inner side of the opalescent screen need cleaning (figs. 18-21).
- Remove the screen by pressing the clips gently with a screwdriver (fig. 20).

Replacing the tube lamps

After approx. 10.000 hours (many years in normal use) we advise you to replace the tube lamps.

Have the tube lamps replaced only by Philips Original Bright Light PL-L 55W.

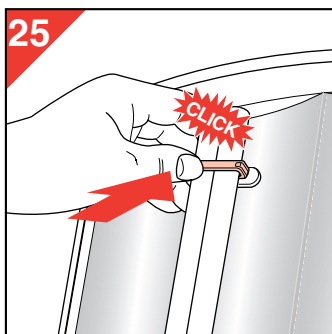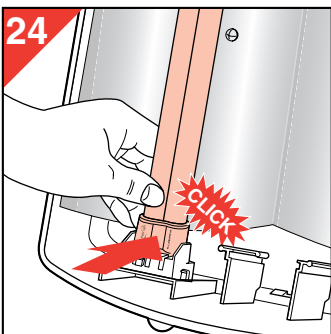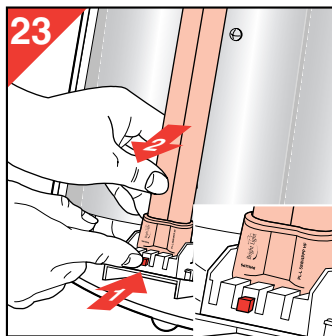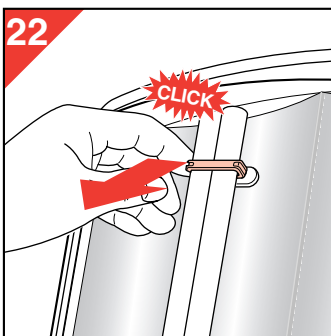You can have the tube lamps replaced by authorized service personnel or replace them yourself. If you wish to replace the lamps yourself:

- Remove the screen (figs. 18 and 19).
 - Bend the bracket, which keeps the lamp to its place, a little sideways and remove it (fig. 22).
 - Press the red button and pull the tube lamp forward, out of the lampholder (fig. 23).
 - When placing a tube lamp, follow these instructions in reverse order (figs. 24 and 25).
- Note:** While sliding the tube light into the lampholder you need not press the red button. As soon as you hear a 'Click!', the lamp is placed correctly.

The tube lamps contain substances which may pollute the environment. When discarding lamps, please ensure that they are separated from the normal household refuse and that they are disposed of at an officially assigned place.


Replacing the mains cord

If the mains cord of this appliance is damaged, it must only be replaced by Philips or their service representative, as special tools and/or parts are required.

This appliance meets the requirements of MDD 93/42/EEC for medical appliances.

Environment

Tube lamps contain substances that may pollute the environment. The lamps of the appliance contain substances that can be harmful to the environment. When you discard the lamps, do not throw them away with the normal household waste but hand them in at an official collection point.

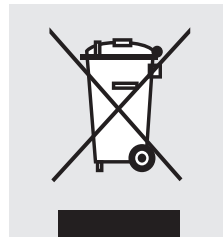- Do not throw the appliance away with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this you will help to preserve the environment.

Guarantee & service

If you need information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

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