# Always here to help you

Register your product and get support at  
[www.philips.com/welcome](http://www.philips.com/welcome)


HP3643

Question?  
Contact  
Philips

# PHILIPS
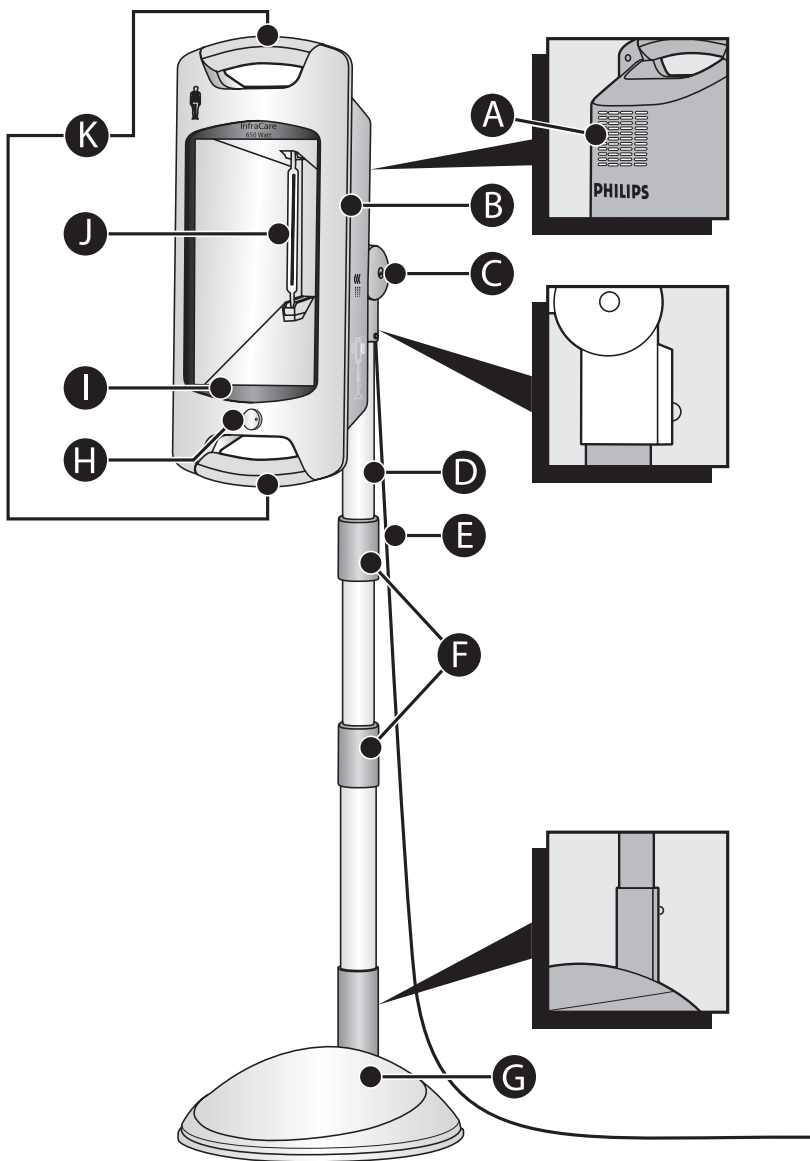
<b>ENGLISH</b>	<b>6</b>
<b>DEUTSCH</b>	<b>14</b>
<b>ESPAÑOL</b>	<b>22</b>
<b>FRANÇAIS</b>	<b>30</b>
<b>ITALIANO</b>	<b>38</b>
<b>NEDERLANDS</b>	<b>46</b>
<b>PORTUGUÊS</b>	<b>54</b>

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

With the InfraCare, you can comfortably treat muscle and joint problems and relieve muscle and joint pain at home. Just like the sun, the appliance emits infrared light. The special filter only passes the type of light required to achieve the intended therapeutic effect. The InfraCare can also be used for relaxation or preparation for a massage.

## General description (Fig. 1)

- A Cooling vents
- B Lamp housing
- C Hinge for tilting lamp housing
- D Pedestal tube
- E Mains cord
- F Height adjustment rings
- G Pedestal
- H On/off switch (I/O) with power-on light
- I Filter
- J Infrared halogen lamp
- K Handgrips

## Important

Read this user manual carefully before you use the appliance and save it for future reference.

### Explanation of symbols

- The symbols on the InfraCare have the following meaning:
  - 1 This symbol warns you against hot surfaces (Fig. 2).
  - 2 This symbol indicates that you have to read the user manual carefully before you use the appliance. Save the user manual for future reference (Fig. 3).
  - 3 This symbol indicates that this appliance is double insulated (Class II) (Fig. 4).
- This appliance meets the requirements of Medical Device Directive MDD93/42/EEC.

### Danger

- Water and electricity are a dangerous combination! Do not use this appliance in wet surroundings (e.g. in the bathroom or near a shower or swimming pool).
- Do not let water run into the appliance or spill water onto the appliance.

### Warning

- Check if the voltage indicated on the appliance corresponds to the local mains voltage before you connect the appliance.
- If the mains cord is damaged, you must have it replaced by Philips, a service centre authorised by Philips or similarly qualified persons in order to avoid a hazard.

- Do not leave the appliance unattended when it is switched on.
- This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- This appliance is not intended for use by children. Children should be supervised to ensure that they do not play with the appliance.
- This appliance is intended for the treatment of muscle and joint pain and muscle and joint stiffness, not for the treatment of abdominal, facial and other sorts of pain.
- To prevent overheating of the skin, avoid very short treatment distances. The warmth should be comfortable throughout the treatment.
- To prevent eye damage, do not look into the light of the lamp when it is on.
- Do not fall asleep during treatment.
- Prevent the treatment area from cooling down too fast immediately after treatment.
- If there is no improvement after 6-8 treatments, stop using the appliance and consult your doctor.

### **Contraindications for InfraCare treatment**

- Do not use the appliance without consulting your doctor first if you suffer from severe diseases such as heart disease, acute inflammatory diseases, thrombosis, blood coagulation disorders, adrenal suppression, systemic lupus erythematosus or malignant diseases.
- Do not use the appliance if you are oversensitive to infrared light or if you have oedema. If in doubt, consult your doctor.
- Do not use the appliance if you are insensitive to infrared radiation and heat due to a disease (e.g. advanced diabetes) or the use of painkillers. If in doubt, consult your doctor.
- Do not use the appliance if you suffer from conditions associated with nerve demyelination such as carpal tunnel syndrome, multiple sclerosis and ulnar neuropathy. If in doubt, consult your doctor.
- Do not use the appliance if you have an implanted device, particularly not if it is located in the vicinity of the treatment area. If in doubt, consult your doctor.
- If you are pregnant, do not use the appliance to treat areas close to your abdomen or your lower back area, to avoid overheating of these areas. If in doubt, consult your doctor.
- Do not use the appliance if you have a swelling or an inflammation, as heat could aggravate the complaints. If in doubt, consult your doctor.

## 8 ENGLISH

- Do not use the appliance on damaged or irritated skin or on wounds. If in doubt, consult your doctor.
- Do not use the appliance to treat areas on which you have recently applied topical counterirritants (creams or ointments containing e.g. menthol), as this can cause a mild inflammatory reaction.

### **Caution**

- If the wall socket used to power the appliance has poor connections, the plug of the appliance becomes hot. Make sure you plug the appliance into a properly connected wall socket.
- Do not subject the appliance to heavy shocks.
- Always unplug the appliance after use and in case of a power failure.
- Place the appliance on a stable, level surface. Make sure there is at least 35cm free space around it to prevent overheating.
- Do not use the appliance at room temperatures lower than 10°C or higher than 35°C.
- Make sure the vents in the back of the appliance remain open during use.
- The appliance is equipped with an automatic protection against overheating. If the appliance is insufficiently cooled (e.g. because the cooling vents are covered), the appliance switches off automatically. Unplug the appliance and remove the cause of overheating. Once the appliance has cooled down, it can be switched on again. To do so, plug the appliance in again and press the on/off switch.
- To avoid getting your fingers stuck when you adjust the height, hold the lamp housing by one of the handgrips or by the hinge as long as the height adjustment rings are unlocked.
- Let the appliance cool down for approx. 15 minutes before you clean or store it.
- Never use the appliance when the filter is damaged, broken or missing.
- If you have just been swimming or have just taken a shower, make sure you dry your skin properly before you use this appliance.

### **Electromagnetic fields (EMF)**

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Using the appliance

### The effects of the InfraCare

The InfraCare produces infrared light, which penetrates deeply into the skin and warms the tissues. The warmth modulates the neuronal signals to the brain and thus reduces pain. Warmth also stimulates the blood circulation and dilates the blood vessels, which accelerates the transport of substances necessary to rebuild and nourish the body's tissues. Furthermore, infrared light speeds up the metabolic process and the removal of waste substances from the body. Because it makes tissue more pliable, warmth also reduces stiffness and makes joints more flexible. Because of these effects, treatment with the InfraCare can temporarily relieve pain due to muscle and joint problems.

*Note: Please read the 'Warning' section in chapter 'Important' and follow the instructions given there before you use the appliance.*

### Applications

#### Therapeutic applications

The InfraCare is a local warmth therapy appliance. Possible therapeutic applications:

- Treatment of occasional or chronic muscle pains, stiff muscles and stiff joints.
- Treatment of lower-back ache.

#### Other applications

You can also use the InfraCare for the following purposes:

- Muscle relaxation
- Preparation for a massage

*Note: For more information about the health benefits, see [www.philips.com](http://www.philips.com), search word InfraCare (Fig. 5).*

### Preparing for use

- 1** Place the pedestal on a stable, level surface (Fig. 6).
- 2** Insert the pedestal tube into the pedestal until it locks into place (Fig. 7).

*Note: There is only one way to attach the tube to the pedestal.*

- 3** Attach the lamp housing to the top of the pedestal tube (Fig. 8).

**To make sure all parts are properly assembled, check if the spring-loaded lock buttons are visible in the holes. (See figure 1).**

- 4** Make sure there is at least 35cm free space around the appliance to prevent overheating (Fig. 9).

Thanks to the pedestal, you can use the appliance in many different positions. Below you find some examples.

- To treat your back from shoulders to hips, place the appliance beside a bed, with the lamp housing in highest position and tilted down. Lie down comfortably on your belly (Fig. 10).
- To treat your back, place the appliance beside a bed, with the lamp housing in lowest position. Lie down comfortably on one side (Fig. 11).
- To treat your upper back and shoulders, place the appliance behind a chair or other seat, with the lamp housing in horizontal position. Sit down comfortably with your back towards the appliance (Fig. 12).
- To treat your back, place the appliance beside a chair or stool, with the lamp housing in vertical position. Sit down comfortably with your back towards the appliance (Fig. 13).

*Tip: Small persons who want to treat their lower back while sitting on a chair can put their feet on a low footstool for more comfort.*

*Tip: Tall persons who want to treat their lower back can lie down on a bed with their back towards the appliance.*

## 10 ENGLISH

- To treat your hip and upper leg, place the appliance beside a chair, with the lamp housing in horizontal position. Sit down comfortably on a chair or stool (Fig. 14).

### **5** Adjust the lamp housing.

**To avoid getting your fingers stuck when you adjust the height, hold the lamp housing by one of the handgrips or by the hinge as long as the height adjustment rings are unlocked.**

*Note: The pedestal tube consists of three parts and two height adjustment rings. You can adjust the height with both rings.*

- To adjust the height of the InfraCare, hold the lamp housing by one of the handgrips or the hinge with one hand while you turn the height adjustment rings with the other hand.
- Turn the height adjustment ring anticlockwise until the 'unlocked' symbol is visible through the hole in the ring (1). Then pull up the lamp housing or push it down to the desired height (2) (Fig. 15).
- When you have set the lamp housing to the desired height, turn the height adjustment ring clockwise until the 'locked' symbol is visible through the hole in the ring. The height adjustment ring is now locked and the height of the lamp housing is fixed (Fig. 16).
- To tilt the lamp housing upwards or downwards, grasp it by its handgrips and move it to the desired position (Fig. 17).
- To put the lamp housing in vertical or horizontal position, grasp it by its handgrips and turn it to the required position (Fig. 18).
- Make sure the pedestal base is horizontal (Fig. 19).

### **6** Sit or lie down at the right distance from the appliance.

- At a distance of 40-50cm to the body, the InfraCare treats an area of about 40cm by 60cm. At a larger distance, the area treated becomes larger and the infrared light and therefore the heat less intense. At a smaller distance, the area treated is smaller and the infrared light and therefore the heat more intense.
- Make sure the distance between your body and the appliance is large enough to prevent accidental contact with the hot filter.
- **Choose a distance at which the heat is comfortable and increase the distance if the heat becomes too intense.**

*Note: The heat of the appliance does not reach its full intensity until several minutes after you switch on the appliance. This does not mean, however, that the appliance is not effective during these first minutes.*

### **Duration of treatment**

---

The duration of the treatment depends on the type of treatment and the person to be treated. Treatment of a single muscle or joint should last approximately 15 minutes (17 minutes in Korea and Taiwan). If necessary, the treatment can be repeated several times a day. A longer treatment time than approximately 15 minutes (17 minutes in Korea and Taiwan) does not lead to better results. You can repeat the course of several treatments a day on several consecutive days to get a result. If there is no improvement after 6-8 treatments, stop using the appliance and consult your doctor. For chronic problems, the appliance can be used for longer periods of time, as long as you follow the instructions given in the 'Warning' section in chapter 'Important'.

### **Using the InfraCare**

---

- 1** Put the mains plug in the wall socket.
- 2** Place the light in the right position and make sure there is enough free space around the light.
- 3** Set the on/off switch to 'I' to switch on the appliance (Fig. 20).
  - ▶ The power-on light goes on (Fig. 21).
  - ▶ The appliance stays on for approx. 15 minutes (17 minutes in Korea and Taiwan) and then automatically switches off.
  - ▶ If you want to use the InfraCare longer, simply press the on/off switch again.

## Cleaning

Unplug the appliance and let it cool down for about 15 minutes before you clean it.

Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.

- 1 Clean the lamp housing, the pedestal and the tube with a moist cloth. Do not let any water run into the appliance.
- 2 You can clean the outside of the filter with a soft cloth sprinkled with some methylated spirits.

## Storage

- 1 Unplug the appliance and let it cool down for about 15 minutes before you store it.
- 2 Turn the height adjustment ring anticlockwise to unlock it (1) and push the lamp housing down to its lowest position (2) (Fig. 15).
- 3 Turn the ring clockwise to fix the lamp housing at this height (Fig. 16).
  - When you lift or carry the appliance, always make sure the lamp housing is in the lowest position. Also make sure that you grab one of the handgrips or the hinge of the lamp housing with one hand and the pedestal tube with the other hand (Fig. 22).

To prevent unintentional unlocking of the height adjustment rings, do not put your hand on a height adjustment ring while you carry or lift the appliance.

- 4 Store the InfraCare in a dry place.

## Replacement

### Filter

The filter protects against exposure to undesired light.

If the filter is damaged or broken, do not use the appliance anymore. Contact the Philips Consumer Care Centre for information.

### Lamp

The effectiveness of the lamp does not diminish in the course of time.

The infrared halogen lamp has a lifetime of approx. 500 hours. If the lamp stops functioning, do not attempt to replace it yourself, but have it replaced by a service centre authorised by Philips. Only they have the knowledge and the skills required for this job and the original spare parts for the appliance.

### Mains cord

If the mains cord is damaged, you must have it replaced by Philips, a service centre authorised by Philips or similarly qualified persons in order to avoid a hazard.

## Environment

- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 23).

## Guarantee and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the separate worldwide guarantee leaflet.

## Troubleshooting

This chapter summarises the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, contact the Consumer Care Centre in your country.

Problem	Possible cause	Solution
The appliance does not go on.	The plug is not inserted into the wall socket properly.	Put the plug in the wall socket properly.
	There is a power failure.	Connect another appliance to check if the power supply works.
	The lamp is defective.	Contact the Philips Consumer Care Centre.
	The mains cord is damaged.	If the mains cord is damaged, you must have it replaced by Philips, a service centre authorised by Philips or similarly qualified persons in order to avoid a hazard.

## Specifications

Model	HP3643
Rated voltage (V) Hong Kong/Singapore/Europe	220-230
Rated voltage (V) China/Korea	220
Rated voltage (V) Malaysia	240
Rated voltage (V) Taiwan	110
Rated frequency (Hz) Hong Kong/Singapore/Europe	50/60
Rated frequency (Hz) China/Malaysia	50
Rated frequency (Hz) Korea/Taiwan	60
Rated power input (W)	650
Insulation class	II
Classification applied part	not classified
Ingress of water	IPX0
IR type	IR A/B/C
Irradiance (W/m <sup>2</sup> )	approx. 1000
Field dimensions (in cm <sup>2</sup> )	60x40
Battery type	not applicable

**Operation conditions**

Temperature	from +10°C to +35°C
Relative humidity	from 30% to 90% (no condensation)

**Storage conditions**

Temperature	from -20°C to +50°C
Relative humidity	from 30% to 90% (no condensation)

Circuit diagrams, component part lists and other technical descriptions are available on request. As user of this appliance, you should take measures to prevent the occurrence of electromagnetic interference or other types of interference between this appliance and other appliances.

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