ENGLISH 4

DEUTSCH 16

FRANÇAIS 29

NEDERLANDS 43

ESPAÑOL 56

ITALIANO 69

PORTUGUÊS 82

TÜRKÇE 94

Introduction

The Cellesse is the result of years of research by the Philips Personal Care Institute. Having investigated the causes of cellulite and tested all the various treatments available, we've now developed the effective way to fight cellulite you're looking for.

This booklet explains what causes cellulite and describes how the Cellesse can help you solve this typically feminine problem.


Important: Before you start using the appliance, please take the time to read this booklet thoroughly (especially the section headed 'Read carefully before using Cellesse').

Cellulite - a typically feminine problem.

- Cellulite is the name given to bumpy, pitted skin that arises when lumpy fat deposits situated just under the skin become visible. Because of its appearance, it's often called "orange peel skin".**

Cellulite only affects certain parts of the body, particularly the hips, thighs and buttocks.

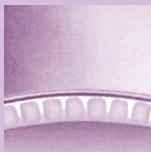- The natural mattress-like structure of women's skin**

Cellulite is the result of two things: the mattress-like structure of women's skin and the tendency to retain fluid in the body.

Men and women have different skin structures. Women's skin naturally has a mattress-like structure. Fat is deposited in 'pockets' that stand at right angles to the surface, giving a structure that resembles a mattress with springs. This mattress-like structure becomes clearly visible when the skin is squeezed. As a result of their skin structure, all women run the risk of developing cellulite.

► Fluid retention

In most cases, the mattress-like structure becomes visible as a result of excess fluid in the fat layer. Women tend to retain more fluid than men, especially in their legs. Excess fluid may accumulate around the fat cells and cause the cellulite layer to swell. If this fluid retention is allowed to continue for months or even years, the fat deposits may thicken and harden. Eventually they'll become hard lumps which can be seen on the surface of the skin even when it's not squeezed.

When it reaches an advanced stage, cellulite may be painful when the skin is squeezed or touched. At that point, cellulite is more difficult to reverse. At earlier stages, however, its effects can be considerably reduced.

Cellulite - what can you do about it?

Once developed, cellulite won't disappear or diminish on its own. On the contrary, it's likely to

increase. The only solution is to make the fight against cellulite part of your regular personal care routine and to begin treatment right away, before the problem gets worse.

The good news is that cellulite can be reduced. It will never disappear completely, however, and will return if you stop doing something about it.

The Cellesse: the cellulite massage system that's been proven to work!

Philips has therefore specially developed the compact and convenient Cellesse to fight cellulite in a effective and comfortable way. The Cellesse gets right down to the cellulite layer within the skin and can easily be used in the privacy of your own home.


Cellesse massages the cellulite layer in your skin in an intensive, pleasant way. This massage mobilises the excess fluid, which is then removed naturally by the body through the normal circulation and lymphatic drainage systems. The result is a reduction in the bumpiness of the cellulite layer and therefore a smoother skin.

The Cellesse in action


The Cellesse has been specially developed for those parts of the body that are prone to cellulite (thighs, hips and buttocks). The appliance works

through a massage system that's been shown in tests to be highly effective - Vacumotion™.

Vacumotion™ works by providing a continuous lifting and rolling movement, which massages the cellulite layer in an intensive and pleasurable way. It consists of two rollers with a suction opening in between.

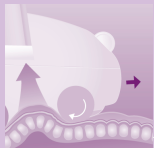- ▶ The front roller gently picks up the skin and guides it into the suction bell.
- ▶ The vacuum lifts the skin even more to form a fold, enabling the massage to reach the cellulite layer effectively.

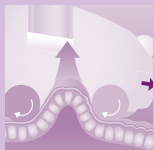- ▶ The two rollers, working together, give the cellulite layer a pleasant and deep massage, removing the main cause of cellulite - the excess fluid surrounding the fat cells. The rear roller then gently guides the skin away.

Following an effective treatment programme for the best results

We recommend you to use the Cellesse three times a week, spending 10 - 15 minutes on each leg each time. Clinical tests have shown that this type of regular treatment can significantly smoothen the cellulite layer in the skin within one month. Naturally, results will vary from individual to individual, depending on how bad the cellulite is to start with and on the individual's lifestyle.

Don't stop the treatment! If you do, the cellulite layer will start to swell again after a while, as a result of the natural tendency of women's bodies to retain fluid. We recommend you to use the appliance three times a week during the first two or three months. After this time you can reduce the frequency and follow a maintenance programme of 15 minutes once a week or 5 minutes three times a week. This should be enough to keep your fluid retention under control.

Proven effective in clinical tests

Although the undesirable consequences of cellulite are visible on the surface, its causes lie deep inside the skin. Independent experts at the University of Brussels carried out two types of clinical tests to establish the effectiveness of the Cellesse.

- One type of test showed that using the Vacuumotion™ system only three times a week smoothens the cellulite layer inside the skin by up to 33% in just one month.

These tests used ultrasound scanning to measure the bumpiness of the cellulite layer - the cause of cellulite.

The results achieved by women using the Cellesse were compared with those of a control group of non-users.

Proven to visibly reduce cellulite

Additional tests showed that the Cellesse is also highly effective on the outside of the skin, producing visibly smoother skin.

Read carefully before using the Cellesse

Do not use the Cellesse on hypersensitive skin, problem skin or areas with skin ailments, irritated skin, varicose veins, scars, spots or wounds. Please consult your doctor before using the Cellesse if you're pregnant or have a heart condition, diabetes or circulation problems.

Only use the appliance to massage the outside of your thighs, hips and buttocks. Don't use it on any other parts of the body!

Side-effects during and after use

After using the Cellesse, your skin may become slightly red and feel rather dry. Don't worry: these are the normal side effects of a massage as the blood circulation improves. They should soon disappear. If necessary, you can apply a body lotion or cream after using the Cellesse to moisturise the skin. If these effects do not disappear within 24 hours or if any other undesirable skin reaction occurs, stop using the appliance and consult your doctor.

If you are prone to bruising or have small red veins clearly visible, you should only use higher settings with care, especially during the first few weeks of use. Whenever necessary, you can

always release the On/Off button to immediately lift the appliance off your skin.

For the best results, the appliance should only be used as described on the following pages.

Important

- ▶ Before using the appliance for the first time, read the section 'Read carefully before using the Cellesse'.
- ▶ This appliance is meant for household use only.
- ▶ Check if the voltage indicated on the mains plug corresponds to the mains voltage in your home before you connect the appliance.
- ▶ Use the appliance only with the cord supplied.
- ▶ Do not use the appliance on a wet skin.
- ▶ The appliance has been specially designed for treatment of thighs, hips and buttocks. Do not use the appliance on the upper arms, the inside of the thighs, the belly or the breasts.
- ▶ Keep the appliance away from water!
- ▶ Keep the appliance away from children.
- ▶ Regular cleaning ensures optimum results and a longer life for the appliance.
- ▶ If the mains cord is damaged, it must be replaced by Philips, a service centre authorised by Philips or similarly qualified persons in order to avoid a hazard.

Using the Cellesse

- 1** Make sure your skin is clean, dry and free of grease or other skincare products.

Do not apply any creams or lotions to the skin prior to the massage.

- 2** Make sure that clothes or long hair do not get caught between the rollers when using the Cellesse. Always tie back long hair before you start using the appliance.

- 3** Put the transformer plug in the wall socket.




- 4** Before starting, set the intensity adjustment knob to the preferred setting.

Higher settings will provide more suction power and therefore intensify the massage.


We advise you to start the first treatment on setting 1. Later, either during the first treatment or when your skin has become used to the appliance, you can gradually increase the intensity. If treatment becomes painful, decrease the intensity. You will achieve the best results when you experience the massage as pleasant. Since everyone's skin is different, you will need to discover which intensity setting suits you best.

Treatment with the Cellesse must, above all, be comfortable.

- 5** Slide your hand under the hand-strap.

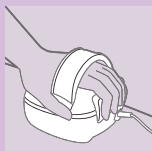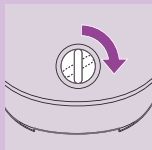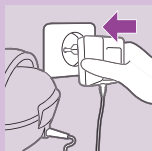
Your fingers should touch the On/Off button. The hand-strap is flexible and adapts itself to your hand.

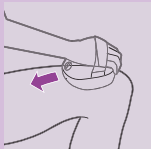
6 Use the Cellesse while you are sitting down, standing up or lying down.

7 Gently position the appliance on your skin, a little above the knee.

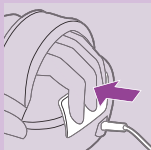8 Press the On/Off button to start the treatment and keep it pressed in.

9 Move the appliance slowly and gradually from a little above the knee to the upper part of your thigh.

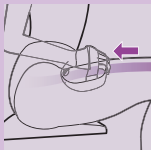
The Vacuumotion massage system will lift the skin, massaging the cellulite layer in an intensive, effective and pleasant way.


10 Always move the appliance in one direction only: upwards towards the heart.

A complete stroke from just above the knee to the upper part of the thigh will take approximately five to seven seconds.

Do not hold the appliance in one place too long
Keep making stroking movements.

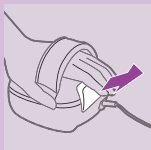11 When you have finished a stroke, release the On/Off button to stop the appliance.

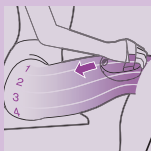
- 12** Lift the appliance off your skin.
- 13** Place the Cellesse back onto your skin, right next to where you began the first stroke, and start a new stroke.
- 14** Repeat this procedure until you have treated your entire thigh, except the inside, with a series of slow, upward, adjacent strokes.

- ▶ If you wish, you may apply a cream, lotion, talcum powder or other type of skincare product to your skin after using the Cellesse.
- ▶ Try to find out which treatment routine suits you best and carry out every treatment with the Cellesse according to this standard routine.
- ▶ The recommended treatment frequency during the first two to three months is three times a week, 10 - 15 minutes on each leg. Read the section headed 'Following an effective treatment programme'.

Cleaning and maintenance

Never immerse the appliance in water or rinse it!

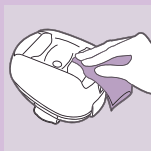
- ▶ For the best results, clean the appliance after every treatment with a moist cloth and, if necessary, a drop of washing-up liquid. Make sure the filter of the suction pump does not get damp.

Do not use alcohol, acetone, petrol, abrasives, scourers, etc., to clean the appliance.


- ▶ The hand-strap can be washed in a washing machine at max. 40°C.

Washing may cause the colour of the strap to change slightly.

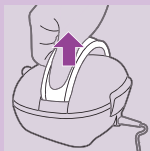- ▶ To remove the hand-strap pull the lip on the inside of the hand-strap and pull the hand-strap upwards at the same time.
- ▶ To fasten the hand-strap again, push the buckles into the recesses and slide them downwards until they snap into position.

Storage

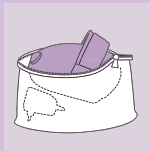- ▶ After use, remove the transformer plug from the wall socket and store the Cellesse in the bag supplied.

Environment

- ▶ Do not throw the appliance away with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this you will help to preserve the environment.

Guarantee & service

If you need information or if you have a problem, please visit the Philips Web site at **www.philips.com** or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

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