PHILIPS

Jordan*

HP 510/610/710/735

English

Teeth for a lifetime with Philips/Jordan

Clear and healthy teeth give the security and confidence of good appearance.

Furthermore, keeping one's teeth for a lifetime is the ultimate goal of dental care. Your own daily dental hygiene as well as the efforts of your dentist and other professionals are aimed at this. Philips/Jordan strongly supports this objective.

What is dental plaque?

To safeguard healthy teeth and gums, proper and frequent removal of plaque in all areas of the mouth is required. Dental plaque is the main cause for dental diseases. It is a soft and sticky film which continually develops on teeth, especially in hard to reach areas as along the teeth line and in between teeth. It has almost the same colour as teeth and is therefore hardly visible. As long as plaque is still soft, it can be easily removed.



What are the effects of plaque?

The bacteria in the plaque produce an acid which can dissolve the enamel of the teeth and the dentine beneath it. This will cause dental decay (caries) (fig. 1)

The bacteria also form substances which may cause inflammation of the gums (gingivitis) or decay of the supporting tissue (periodontitis).

Inflamed gums are generally red and swollen and often bleed when the teeth are brushed (fig. 2).

Thorough and regular removal of plaque in all areas of the mouth will prevent gingivitis from starting and will help stop the spread of periodontitis.

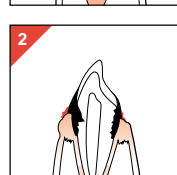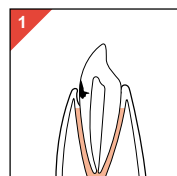Philips/Jordan 2-Action plaque remover

2-Action Brush

The power of this new developed plaque removal system is in the double performance of the oscillating brush. Fast oscillating movements combined with the additional active tip

movements in the space between the teeth, result in a thorough and highly effective cleaning of every area of the teeth (fig. 3).

Controlled Pressure

Proper cleaning technique requires excellent handling of a brush as it is a balancing act between applied pressure and movements.



Too much pressure may damage both teeth and gums, causes unnecessary wear of the brush and also reduces the effectiveness of the brushing operation.

The Controlled Pressure system helps you to control applying the correct pressure when brushing your teeth.

When the pressure level is exceeded, the brush actively warns you by a sudden change in the position of its neck.

This may occur often when first time used.

Clinical tests show that the Controlled Pressure will assist you in developing effective electrical brushing techniques for cleaning at the right pressure, within a few days (fig. 4).

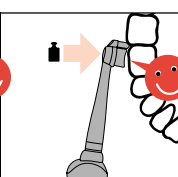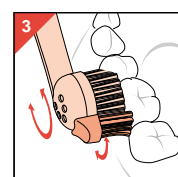Important

- The appliance must always be mounted to the wall.
- First charging: at least 16 hours.
- The mains cord of this appliance cannot be replaced.
- If the cord is damaged, the appliance should be discarded.
- For optimum brushing efficiency and for hygiene: replace the 2-Action brush heads as soon as the hairs begin to stand out or at least every three months.
- It is advisable to use your electrical toothbrush twice a day for at least two minutes.
- The built-in rechargeable battery contains substances which may pollute the environment.
- When eventually discarding your appliance please remove the built-in rechargeable battery (fig. 22-24) and dispose of it at an officially assigned refuse site.
- Please refer to section "Protecting the environment".

General Description (fig. 5)

- 2-Action brush head
- Oscillating brush
- Active tip
- Coloured identification ring (click-on)
- Drive shaft
- "Controlled Pressure" bending head
- Brushing time indicator (Type HP 735 only)
- Speed selector (button)
- On/off switch (button)
- Charge control (Types HP 710 and HP 735 only)
- Screws and rawlplugs for wall mounting
- Charge unit with mains cord
- Brushes storage compartment - for two 2-Action brush heads (Type HP 735 only)
- Inner container for
- Brush peg - for one 2-Action brush head (Types HP 510, HP 610 and HP 710 only)
- Cord storage compartment

Bathroom use

- Using the mounting plate, mount the charge unit to a wall where it cannot accidentally fall into water.
- As the toothbrush itself is waterproof you may choose to use it while having a shower or bath.

Cord storage (fig. 6)

Wall mounting (fig. 7-8a)

How to charge (fig. 9-11)

- Before first use charge for at least 16 hours.
- In normal use, the toothbrush is recharged after each use by just placing it back into the connected charge unit.
- As long as the battery is being (re)charged the green pilot light will be on.

100% recycled paper

4203 000 42725

To guarantee a long life, it is important to let the toothbrush run until it is completely empty and then recharge it for 12 hours at least once every six months.

How to use (fig. 12-16)

- Let the toothbrush do the work for you: just guide the bristle carefully along the teeth. The active tip will find its own way and enter into the spaces in between the teeth. Ensure that the brush touches the entire tooth surface as well as a small section of the gums.
- If the head clicks backwards during brushing, this means that you have applied too much pressure (fig. 17).
- Also refer to the provided Tooth Brushing Instructions card.
- Position (push and turn) the 2-Action brush head on to the appliance (fig. 12).
- Press button (I) to switch on (fig. 14).
- The toothbrush will start at high speed.
- To change speed press button (H) (fig. 15). (Types HP 710 and HP 735 only.)

Minimum brushing time: 2 minutes (fig. 18)

Two minutes is regarded to be the minimum time required for effective plaque removal and gum massage. For even better results you may decide to brush longer. Remember that fluoride in toothpaste also needs some time to act. (Type HP 735 only.) To help you brush systematically (fig. 18) the brushing time indicator (G) will come on for 10 seconds at intervals of 30 seconds.

After brushing your teeth (fig. 19-21)

Cleaning

- Switch the toothbrush on and rinse it clean in lukewarm (max. 60 °C) water.
- The charge unit can be cleaned with a damp cloth. Unplug before cleaning!

(Type HP 735 only.) The inner container is removable and can be washed (also dishwasher safe).

Replacement 2-Action brush heads

For replacement, new 2-Action brush heads are available from your dealer under type n°. HP 790.

Protecting the environment

To remove the built-in rechargeable battery when eventually discarding your appliance: see fig. 22-24.

- Please note: You can also take the appliance to a Philips Service Centre. The staff will be happy to help you.
- Before removing the battery should be completely empty of charge. Let the appliance run until the motor stops.
- Please note: Only unscrew the bottom for battery removal. The toothbrush may not be watertight anymore after the bottom has been unscrewed. For that reason the guarantee is void after opening the toothbrush.
- An undamaged sticker at the bottom side indicates that the appliance has not been opened.

For U.K. only:

Mains plug

The mains plug of the charging unit meets the latest European standards. The plug is completely safe, although the standard UK shaver socket holes may be slightly too wide for the pins of the plug. Please make sure a proper connection is made.

If the plug must be inserted into a 3-pin socket, use a suitable shaver adapter.

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