HX 1525/20




ENGLISH	4
FRANÇAIS	11
DEUTSCH	18
NEDERLANDS	25
ITALIANO	32
ESPAÑOL	39
PORTUGUÊS	46
TÜRKÇE	53

Important

- ▶ Be sure to check if the wall socket to which you connect the charging unit is live. In some bathroom cabinets, the electricity supply to the built-in socket may be interrupted as soon as the light is switched off.
- ▶ Charge the appliance for at least 16 hours before you use it for the first time.
- ▶ The mains cord of the charging unit cannot be replaced. If the cord is damaged, the unit should be discarded. Please apply to your Philips dealer or to the Philips organisation in your country (refer to the World Wide Guarantee leaflet).
- ▶ The built-in battery contains substances that can be harmful to the environment. If you discard the appliance in due course, remove the built-in battery and hand it in at an official collection point or dispose of it in accordance with the regulations. See section "Environmental aspects".
- ▶ ***If your appliance is marked with  the following text applies:***
To protect the environment, this appliance has been equipped with a nickel-hydride (NiMH) battery which contains less than 0.01% cadmium. Still we advise you to make sure sure that, when you discard the appliance, you keep the battery separate from normal household refuse and dispose of it at an officially assigned collection point. Please refer to section "Environmental aspects".
- ▶ ***If your appliance is marked with  the following text applies:***
The built-in rechargeable (NiCd) battery contains substances which may pollute the environment. When you discard your appliance in due course, please ensure that the battery is kept separate

from normal household refuse and disposed of at an officially assigned collection point. Please refer to section "Environmental aspects".

Environmental aspects

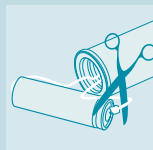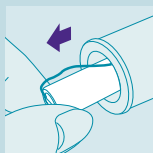
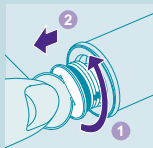
If you discard the appliance in due course, remove the battery and dispose of it in an environmentally safe way.

- ▶ NB: You can also hand the appliance in at a Philips Service Centre, which will dispose of the battery for you.
- ▶ Before you remove the battery, it must be completely empty. Let the appliance run until it stops before removing the battery.
- ▶ NB: Only unscrew the bottom of the appliance if you want to remove the battery.

After the bottom has been unscrewed, the appliance may not be watertight anymore. Therefore the guarantee will become invalid if the bottom of the appliance has been unscrewed. If undamaged, the sticker on the bottom of the appliance will prove that the appliance has not been opened.

Directions for use


Mounting the appliance on the wall

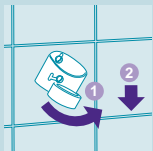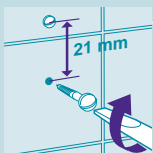1 Mount the charging unit to the wall.

- Do not allow the charging unit to come into contact with water. The toothbrush, however, is water tight and can therefore be used safely anywhere.

Charging

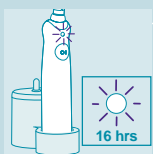
1 Leave the appliance to charge for at least 16 hours before you use it for the first time.


- The toothbrush can be put back in the charging unit after use.
- Only type HX 1525:* If the battery in the toothbrush is being charged, the pilot light is green.

Tip: The toothbrush can be left to charge continuously. When you start using the appliance, we recommend you to use it until the battery is empty and to recharge it completely. Do so three times in a row. You can do this by using the appliance in the normal way without recharging it between brushing sessions until the motor almost stops running. To extend the lifetime of the battery, repeat this procedure every six months.


How to use

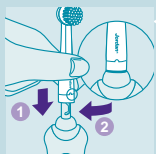- 1 Mount the brush head onto the toothbrush by inserting it into the bayonet catch and turning it. This will make it impossible for the brush head to become detached due to vibration.

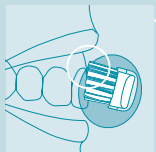- 2 Place the brush in your mouth.

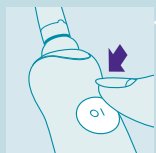
- 3 Press the on/off button to switch the toothbrush on.

Brushing


- 1 Place the brush against the tooth at an angle, in such a way that it covers the surface of the tooth as well as the gum line. Hold the brush in that position for a few seconds and then move on to the next tooth.

- ▶ The round brush head removes plaque from the teeth and beneath the gum margin.
- ▶ The individual moving Active Tip cleans between the teeth.

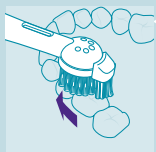
- 2 The brush will do the work for you. All you have to do is guide the brush carefully along the teeth and molars.

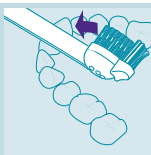
- ▶ Make sure that you reach all surfaces of the teeth (inside, outside and chewing surface) as well as the gums properly. A fixed brushing pattern will prevent you from skipping certain areas.
- ▶ See also the brushing instructions coming with the appliance and the section 'taking optimal care of your teeth'.

Safe brushing

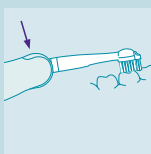- 1 The Controlled Pressure system will keep a check on the pressure you exert on the brush.

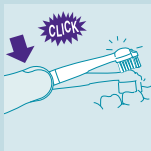
- 2 If you press the brush too hard against your gums, you will feel and hear a click when the brush head flexes backwards. As soon as you reduce the pressure, the brush will return to its normal position.
- ▶ This will probably happen regularly when you first start using the brush, but after some time you will have discovered the optimal brushing pressure, and the brush will flex backwards less often.

Tip: If plaque is removed properly, the gums will have a pink colour and firmly surround the teeth. Healthy gums do not bleed when you eat or brush your teeth.


If, nevertheless, your gums bleed slightly the first few times you use the brush, this bleeding will generally stop after a few days. If it does not, consult your dentist or dental hygienist.

Minimum brushing time: 2 minutes

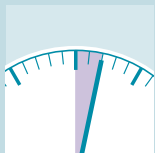
- ▶ In general, 2 minutes will be enough to clean your teeth properly.
- ▶ Naturally, you can also brush longer than 2 minutes.


After brushing

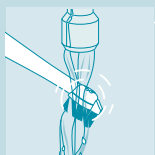
- 1 Let the toothbrush run while you rinse the brush head under the tap. Use lukewarm water (max. 60°C).

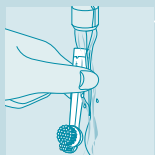- 2 Also rinse the inside of the brush head.

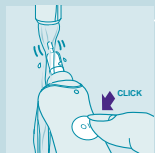
- 3 Switch the toothbrush on again to rinse toothpaste residues from the drive shaft.

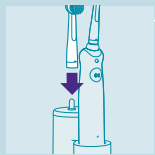- 4 Store the brush head by placing it onto the charging unit.
- 5 You can clean the charging unit with a moist cloth.
 - ▶ First unplug the charging unit!

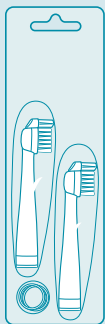Brush head replacement

- 1** Replace the brush head every three months for the best brushing results.
 - ▶ The brush head should also be replaced when the bristles start to bend outwards.
 - ▶ New brush heads are available from the shop in which you have purchased your Philips Jordan toothbrush.

Tip: Worn bristles will move plaque around rather than remove it. They may also damage the gums. It is therefore important that you replace the brush head every three months.

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