
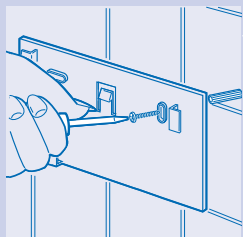
ENGLISH	4
FRANÇAIS	17
DEUTSCH	31
NEDERLANDS	45
ITALIANO	59
ESPAÑOL	73
PORTUGUÊS	86
TÜRKÇE	101

Mounting the appliance on the wall



1 Mount the charging unit on the wall by means of the mounting plate. Choose a place where the appliance cannot fall into the water accidentally.

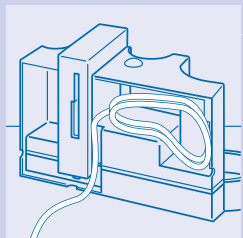- ▶ The toothbrush itself is waterproof, so you may use it while having a shower or bath.

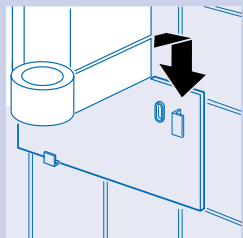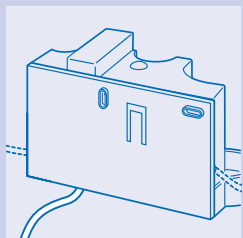2 Push excess mains cord into the storage space in the back of the charging unit.

Important

- ▶ Charge the toothbrush for at least 16 hours before you use it for the first time.
- ▶ The mains cord of the charging unit cannot be replaced. If the cord is damaged, the unit should be discarded. Turn to your Philips dealer or to the Philips organisation in your country for help or advice (refer to the World Wide Guarantee leaflet).
- ▶ The mains plug of the charging unit meets the latest European standards. The plug is safe, although the standard UK shaver socket holes may be slightly too wide for the pins of the plug. Make sure a proper connection is made.
- ▶ Use a suitable adapter if you have to connect the appliance to a standard UK 3-pin socket.
- ▶ If your toothbrush is marked with **(NiMH)** the following text applies:
To protect the environment, this toothbrush has been equipped with a nickel-hydride (NiMH) battery which contains less than 0.01% cadmium. Nevertheless we advise you to make sure that, when you discard the toothbrush, you keep the battery separate from normal household refuse and dispose of it at an officially assigned collection point. Please refer to the section "Environmental aspects".
- ▶ If your toothbrush is marked with **(NiCd)** the following text applies:
The built-in rechargeable (NiCd) battery contains substances which may pollute the environment.
When you discard your toothbrush in due course, please make sure that the battery is kept separate from normal household refuse and dispose of it at an officially assigned collection point. Please refer to the section "Environmental aspects".
- ▶ When you use the Oral Water Jet for the first time or when you have not been using it for a long time, it is necessary to set the volume control to its maximum setting (5) and the jet tip nozzle to position • (for single-jet operation) to start the water flowing to the jet tip. If the water doesn't flow out after 10 seconds, put the jet tip nozzle in your mouth and suck it shortly.

Taking optimal care of your teeth

Healthy teeth are to be treasured, which is logical, for we want them to last a lifetime. As plaque is the major cause of tooth decay, proper dental care is essential. By brushing carefully every day you will remove plaque and keep your teeth and gums healthy longer

What is plaque?

Plaque is a soft, sticky layer on our teeth. It forms mainly on hard-to-reach areas such as interdental spaces and along the gum line. Plaque mainly consists of bacteria that turn the sugars in our food into harmful acids. Due to its yellowish white colour, plaque is hard to detect on the teeth. Fortunately, as long as the plaque is still soft, it can easily be removed by means of careful and correct brushing.

What are the consequences of plaque?

The acids in plaque dissolve the enamel on our teeth and molars, which may lead to the development of cavities. The bacteria produce toxins which may cause the gums to become inflamed (gingivitis). Inflamed gums are often red

and swollen and may start bleeding during brushing or eating. Moreover, the inflammation may spread from the gum line to the jaw bone beneath, which may eventually cause the teeth and molars to fall out. If plaque is not removed daily, it may calcify into tartar adhering strongly to the teeth. Tartar, in its turn, provides an excellent platform for new plaque growth, thus further increasing the risk of caries and inflammation.

What is the best toothbrushing method?

Toothbrushing is an activity that needs to be carried out with care. Therefore, take your time to reach all areas. Electric toothbrushes have been proven to be more effective than manual toothbrushes when it comes to removing plaque. With the Philips Jordan Sensiflex, brushing each side of the tooth for a few seconds is enough to clean both the tooth and the surrounding gums. In this way, you can clean all your teeth in a few minutes (see also the brushing instructions). If you like, you can ask your dentist or dental hygienist for a personal brushing advice.

Correct positioning of the toothbrush

The Philips Jordan Sensiflex has been designed to make the type of movement that best removes plaque. As there is no need for you to move the brush to and fro, you can focus on placing the brush on your teeth correctly. The quickly oscillating bristles of the brush first detach and then wipe away the plaque. The extra long and soft bristles of the Active Tip ensure thorough and effective cleaning between the teeth and at the back of the mouth.

Enough time

By guiding the brush slowly along the insides, outsides and chewing surfaces of all teeth and molars, you will clean all areas in which plaque is likely to develop. Hard-to-reach areas such as behind the molars, along the gum line and behind the front teeth, require extra attention. It is therefore important that you take enough time for your daily oral hygiene.

Optimal pressure

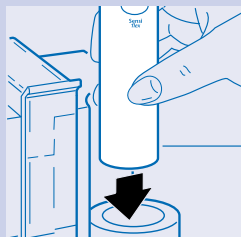
Brushing too hard may cause damage to your teeth and gums. And, contrary to what you might expect, it removes plaque less effectively than gentle brushing. The unique Controlled Pressure system will help you keep a check on your brushing pressure. If too much pressure is applied to the brush, the toothbrush will warn you immediately by means of a click that you will both hear and feel. Then the brush head will flex back. Clinical tests have shown that, thanks to this Controlled Pressure system, brushing with the Philips Jordan Sensiflex results in a better and safer brushing technique.

Oral Water Jet

The Oral Water Jet is an additional device to assist you in the daily care for your teeth. It gives your mouth a delightfully fresh feeling, particularly when mouthwash has been added to the water.

Directions for use: toothbrush


Charging the toothbrush


- 1 Leave the toothbrush to charge for at least 16 hours before you use it for the first time.
 - ▶ The toothbrush can be put back in the charging unit after use.
 - ▶ If the battery in the toothbrush is being charged, the pilot light is green (only type HX2740).
 - ▶ If the battery is almost empty, the pilot light is red and blinks for 10 seconds after the toothbrush has been switched off (only type HX2740).

Tip: The toothbrush can be left to charge continuously. However, to extend the life of the battery, it is recommended to let the toothbrush run until it is completely empty and then fully recharge it every six months.

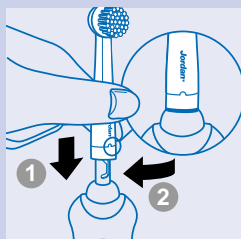
How to use the toothbrush

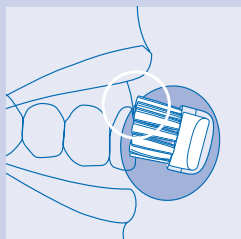
- 1 Mount the brush head onto the toothbrush by inserting it into the bayonet catch and turning it. This will make it impossible for the brush head to become detached due to vibration.

- 2 Place the brush in your mouth.

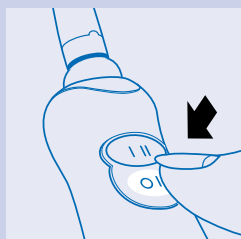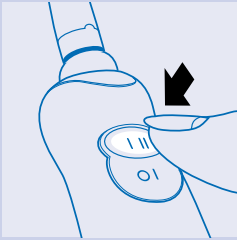
- 3 Press the on/off button to switch the toothbrush on.


- 4** The toothbrush will start oscillating at the highest speed. To switch to a lower speed, press the 2-speed button (only type HX2740).

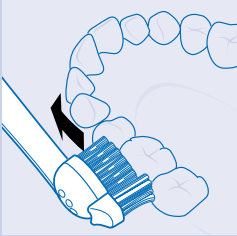
Brushing

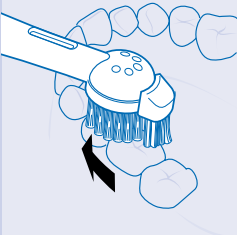- 1** Place the brush against the tooth at an angle, in such a way that it covers the surface of the tooth as well as the gum line. Hold the brush in that position for a few seconds and then move on to the next tooth.
- ▶ The round brush head removes plaque from the teeth and beneath the gum line.
 - ▶ The individually moving Active Tip cleans between the teeth.

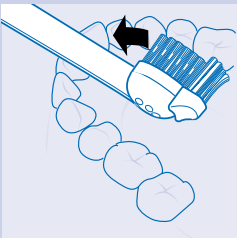
- 2** The brush will do the work for you. All you have to do is guide the brush carefully along the teeth and molars.
- ▶ Make sure that you reach all surfaces of the teeth (inside, outside and chewing surface) as well as the gums properly. A fixed brushing pattern will prevent you from skipping certain areas.


- ▶ See also the brushing instructions coming with the appliance and the section 'taking optimal care of your teeth'.

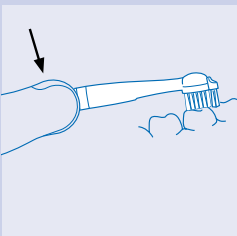Safe brushing

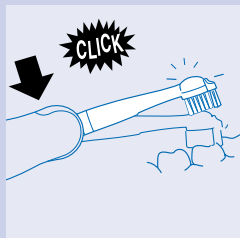- 1** The Controlled Pressure system will keep a check on the pressure you exert on the brush.


- 2 If you press the brush too hard against your gums, you will feel and hear a click when the brush head flexes backwards. As soon as you reduce the pressure, the brush will return to its normal position.
 - ▶ This will probably happen regularly when you first start using the brush, but after some time you will have discovered the optimal brushing pressure, and the brush will flex backwards less often.

Tip: If plaque is removed properly, the gums will have a pink colour and firmly surround the teeth. Healthy gums do not bleed when you eat or brush your teeth (see also the section 'taking optimal care of your teeth').


If, nevertheless, your gums bleed slightly the first few times you use the brush, this bleeding will generally stop after a few days. If it does not, consult your dentist or dental hygienist.

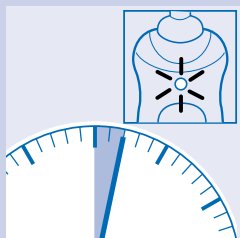Minimum brushing time: 2 minutes


- ▶ In general, 2 minutes will be enough to clean your teeth properly.
- ▶ Therefore the brush is equipped with a red pilot light that will blink for 40 seconds when you have been brushing for 2 minutes (only type HX2740).
- ▶ Naturally, you can also brush longer than 2 minutes.

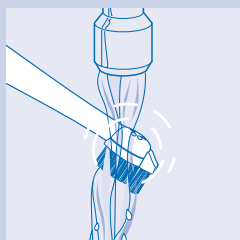
After using the toothbrush

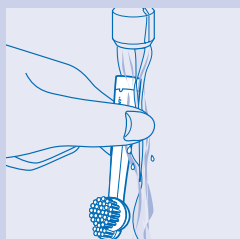
- 1 Let the toothbrush run while you rinse the brush head under the tap. Use lukewarm water (max. 60°C).

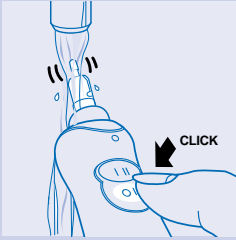- 2 Also rinse the inside of the brush head.

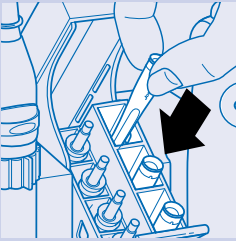
- 3 Switch the toothbrush on again to rinse toothpaste residues from the drive shaft.

- 4 Store the brush head upside down in the storage compartment, with the coloured ring on top.
 - ▶ You can clean the charging unit with a moist cloth.
 - ▶ First unplug the charging unit!

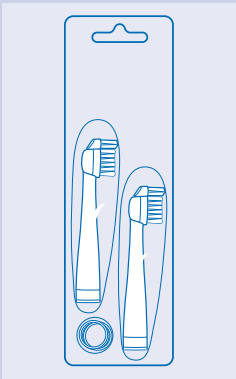
Brush head replacement

- 1 Replace the brush head every three months for the best brushing results.
 - ▶ The brush head should also be replaced when the bristles start to splay.
 - ▶ New brush heads are available from the shop in which you have purchased your Philips Jordan toothbrush.


Tip: Worn bristles will move plaque around rather than remove it. They may also damage the gums. It is therefore important that you replace the brush head every three months.

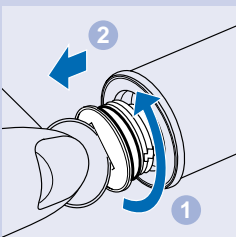Environmental aspects


When you discard the toothbrush in due course, remove the battery and dispose of it in an environmentally safe way.

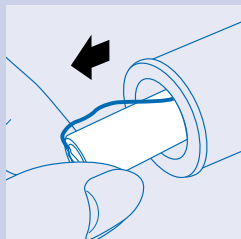- ▶ NB: You can also hand the toothbrush in at a Philips Service Centre, which will dispose of the battery for you.
- ▶ Before you remove the battery, it must be completely empty. Let the toothbrush run until it stops before removing the battery.


► NB: Only unscrew the bottom of the toothbrush if you want to remove the battery.

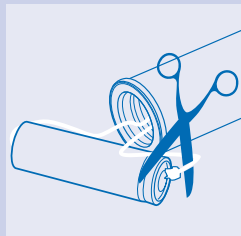After the bottom has been unscrewed, the toothbrush may not be watertight anymore. Therefore the guarantee will become invalid if the bottom of the toothbrush has been unscrewed. If undamaged, the sticker on the bottom of the toothbrush will prove that the toothbrush has not been opened.


Directions for use: Oral Water Jet

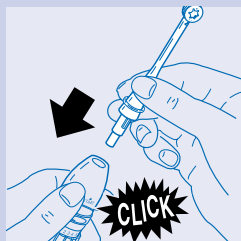
When you use the Oral Water Jet for the first time or when you have not been using it for a couple of days:


1 Click the jet tip onto the appliance.

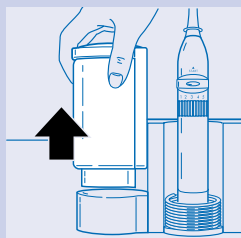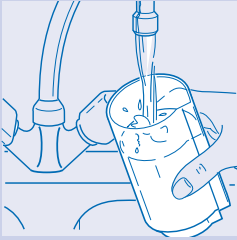
2 Remove the water reservoir from its base.


3 Fill the water reservoir with water.

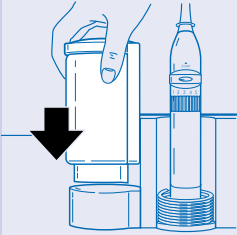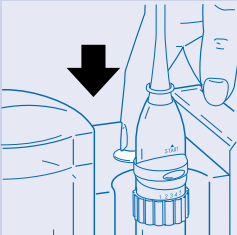4 Put the water reservoir back onto its base. Press it home properly.

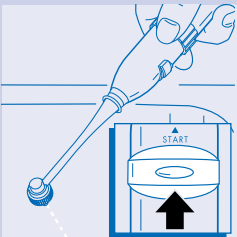5 Press the power button to switch the Oral Water Jet on.


6 Push the switch slide on the handgrip upwards to start the water running.

- ▶ When you use the Oral Water Jet for the first time, let the appliance run until the water reservoir is completely empty. Refill the water reservoir.

Tip: You may add mouthwash for extra freshness to the water. After using mouthwash, fill the water reservoir with clean water and pump all the water out of the appliance and into the washbasin.

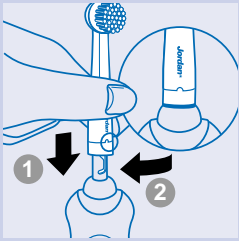
How to use the Oral Water Jet


1 Select position • for cleaning.

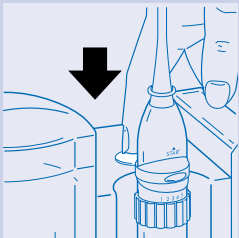
2 Bend over the washbasin and hold the jet tip close to your teeth.

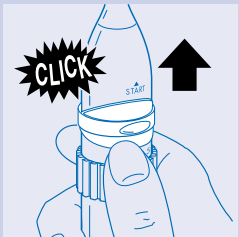3 Press the power button to switch the Oral Water Jet on.

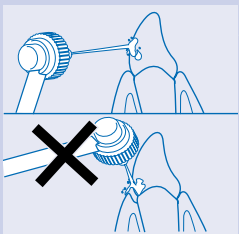
4 Push the switch slide upwards to start the water running.


5 Guide the water jet along the front and the back sides of the teeth and along the gums, and direct it especially at the spaces between the teeth to remove any remaining food particles.

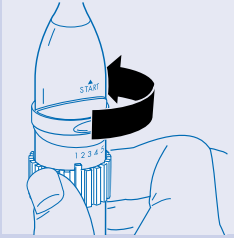- ▶ Make sure that the water jet is at right angles to the teeth.
- ▶ Do not direct the water jet at the gums in such way that the jet penetrates beneath the gum line.
- ▶ You can switch the water supply off with the switch slide before changing the position of the jet tip in your mouth.

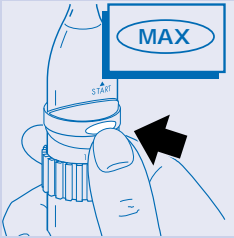
- 6 Adjust the water volume by turning the control ring.


- 7 Press the MAX button if you want to use maximum jet power for a short time (especially when cleaning under bridgework).
 ▶ Release the button to switch back to the set water volume.

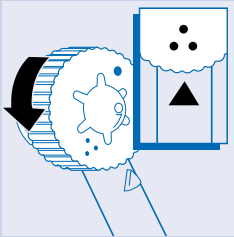Tip: We advise you to use a low water volume setting while you are getting accustomed to the Oral Water Jet and when you are troubled by gum inflammation.

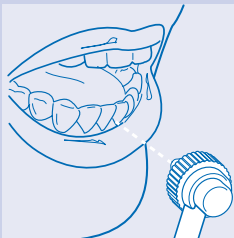
- 8 After cleaning, set the nozzle to position ●● for gum massage.


- 9 Guide the water over your teeth and gums.

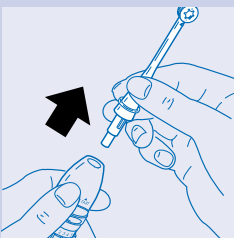After using the Oral Water Jet

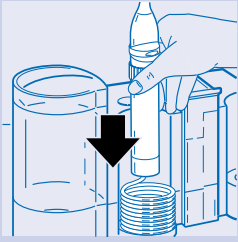- 1 Remove the jet tip and put it in the storage compartment.


2 Put the handgrip back into its holder.

Tip: You can obtain extra or replacement jet tips at your nearest Philips Service Centre.


Cleaning the Oral Water Jet

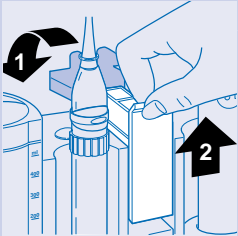- ▶ Remove the mains plug from the wall socket.
- ▶ Clean the appliance with a damp cloth and with some detergent, if necessary.
- ▶ Do not use abrasive or aggressive cleaners, scourers, alcohol, acetone, etc.


1 Clean the storage compartment, the handgrip holder and the water reservoir in hot soapy water.

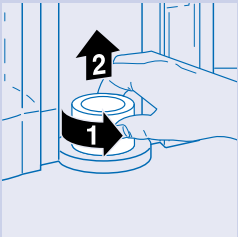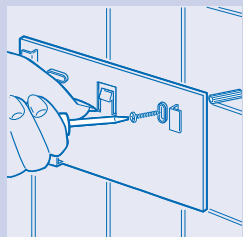- ▶ Rinse them under a running tap afterwards.
- ▶ These items can also be cleaned in a dishwasher

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