# PHILIPS

sensiflex  
dental center HX2225

**HX2225****ENGLISH 4****DEUTSCH 10****FRANÇAIS 17****NEDERLANDS 24****ESPAÑOL 30****ITALIANO 37****PORTUGUÊS 43**

## Important

- ▶ If the mains cord of this appliance is damaged, it must always be replaced by Philips or a service centre authorised by Philips, in order to avoid hazardous situations.
- ▶ The mains plug of the charging unit meets the latest European standards. The plug is safe, although the standard UK shaver socket holes may be slightly too wide for the pins of the plug. Make sure a proper connection is made.
- ▶ When you use the Massage Jet for the first time or when you have not been using it for a long time, it is necessary to set the volume control to its maximum setting to start the water flowing to the jet tip. If the water doesn't flow out after 10 seconds, put the jet tip nozzle in your mouth and suck it shortly.


## Taking optimal care of your teeth

Healthy teeth are to be treasured, which is logical, for we want them to last a lifetime. Healthy gums play an important role. By massaging with the Massage Jet you keep your gums healthy.

Proper dental care is essential. By brushing only, you may not remove all remaining food particles. This may cause caries. By using the Massage Jet, you are able to remove all remaining food particles. Even hard to reach places, such as between your teeth, can be cleaned.

## Preparing the Massage Jet for use

### Mounting the base unit to the wall

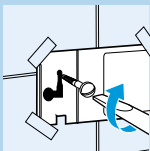
- 1** Mount the base unit to the wall, using 2 screws (3.5x30mm) and 2 wall plugs (S5) (not included).

To drill the holes in the right places, you can use the drilling stencil on the last page of these directions for use.

Do not allow the base unit to come into contact with water.


- 2** Push excess mains cord into the storage space in the back of the base unit.


### Using the Massage Jet

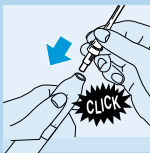
When you use the Massage Jet for the first time or when you have not been using it for a couple of days:

- 1** Click the jet tip onto the appliance.
- 2** Remove the water reservoir from the base unit.

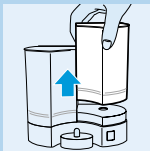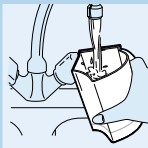
**3** Fill the water reservoir with water.

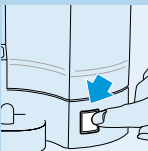**4** Put the water reservoir back onto the base unit. Press it home properly.

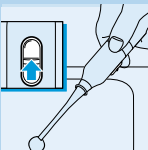
**5** Press the power button to switch the Massage Jet on.

**6** Push the slide on the handgrip upwards to make the water start running.

When you use the Massage Jet for the first time, let the appliance run until the water reservoir is completely empty. Refill the water reservoir.


Tip: You may add mouthwash for extra freshness to the water. After using mouthwash, fill the water reservoir with clean, unmixed water and pump all the water out of the appliance and into the washbasin.

## How to use the Massage Jet

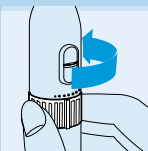
**1** You can adjust the water volume by turning the control ring.


Tip: we advise you to use a low water volume setting while you are getting accustomed to the oral water jet and when you are troubled by gum inflammation.

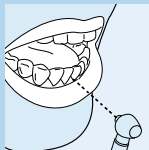
- 2** Bend over the washbasin and hold the jet tip close to your teeth.

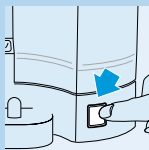- 3** Press the power button to switch the Massage Jet on.

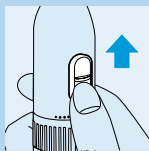
- 4** Push the switch slide upwards to start the water running.

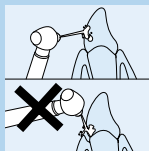
- 5** Guide the water jet along the front and the back sides of the teeth and along the gums, and direct it especially at the spaces between the teeth to remove any remaining food particles.

- ▶ *Make sure the water jet is at straight angles to the teeth.*
- ▶ *Do not direct the water jet at the gums in such way that the jet penetrates beneath the gum line.*

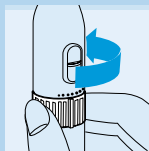- ▶ You can switch the water supply off with the switch slide before changing the position of the jet tip in your mouth.

- 6** Adjust the water volume by turning the control ring.

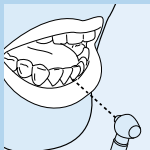
- ▶ **Tip:** We advise you to use a low water volume setting while you are getting accustomed to the Massage Jet and when you are troubled by gum inflammation.


- 7** After cleaning, lower the water volume for gum massage.
- 8** Guide the water over your teeth and gums.

### After using the Massage Jet

---

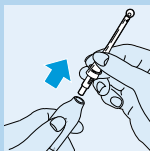- 1** Remove the jet tip and put it in the storage compartment.

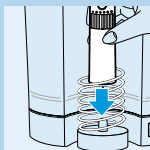- 2** Put the handgrip back into its holder.

**Tip:** You can obtain extra or replacement jet tips at your nearest Philips Service Centre.

### Cleaning the Massage Jet

---

- 1** Remove the mains plug from the wall socket.
- 2** Clean the appliance with a damp cloth (if necessary with some detergent).

Do not use abrasive or aggressive cleaners  
scourers, alcohol, acetone, etc.

- 3** Clean the storage compartment, the handgrip holder, the insert container and the water reservoir in hot soapy water.
  - 4** Rinse them under a running tap.
- The storage compartment, the handgrip holder and the water reservoir may also be cleaned in a dishwasher.

## Guarantee & service

If you need information or if you have a problem, please visit the Philips Web site at [www.philips.com](http://www.philips.com) or contact the Philips Customer Care Centre in your country (you will find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, turn to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

Extra tips about the Massage Jet (for replacement or extension) are available as an option at Philips Service Centres.

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