
# Sonicare For Kids rechargeable sonic toothbrush

**PHILIPS**  
**sonicare**

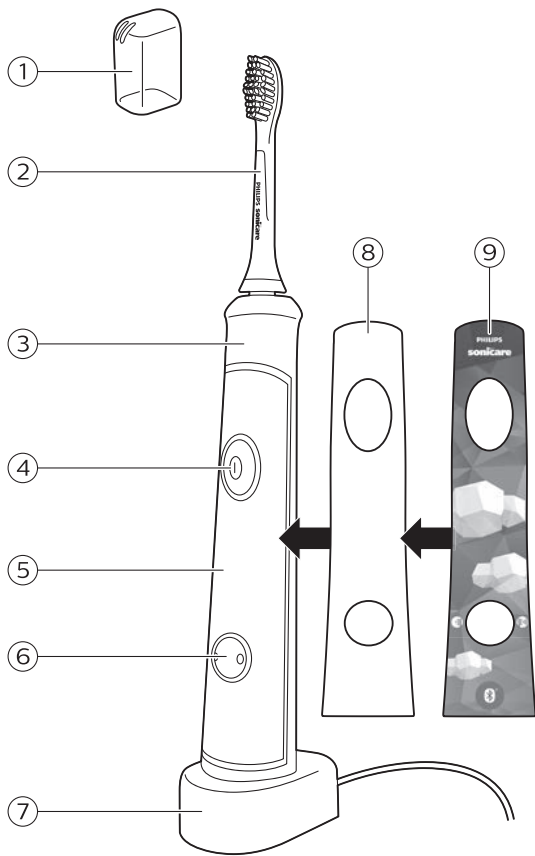English	6
Български	19
Čeština	32
Eesti	45
Hrvatski	58
Magyar	71
Қазақша	85
Lietuviškai	98
Latviešu	111
Polski	123
Română	137
Русский	150
Slovensky	164
Slovenščina	177
Srpski	189

## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at [www.philips.com/welcome](http://www.philips.com/welcome).

You and your child now own a Sonicare connected toothbrush. The Bluetooth connected toothbrush interacts with a fun brushing mobile device application to keep your child engaged while learning proper brushing techniques.

## General description (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head
- 3 Handle with soft grip
- 4 Power on/off button with battery charge indicator
- 5 Removable panel with replaceable stickers
- 6 Brushing mode button with mode lights
- 7 Charger
- 8 Panel
- 9 Sticker

## Important safety information

Read this important information carefully before you use the device and save it for future reference.

### Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water

or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

- This appliance is not a toy. Check brush heads regularly for cracks. Cracked brush heads may present a choking hazard.

## **Warning**

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Always have the charger replaced with one of the original type in order to avoid a hazard.
- Do not use the charger outdoors or near heated surfaces.
- If the appliance is damaged in any way (brush head, toothbrush handle or charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee and support').
- This appliance can be used by children aged from 8 years and above and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they

understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.

## Caution

- Do not clean the brush head, the handle or the charger in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use the Philips Sonicare.
- The Philips Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Philips Sonicare.

- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- The Philips Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.
- Do not charge your phone or tablet in the bathroom.

### **Electromagnetic fields (EMF)**

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Preparing for use

### Applying the replaceable panel

The Sonicare For Kids comes with 1 replaceable panel.

- 1 Remove the adhesive backing from the panel.
- 2 Apply the panel to the white section on the front of the handle (Fig. 2).

Note: Make sure you line up the cut-out areas on the panel with the buttons on the front of the handle.

- 3 Apply your favorite sticker to the panel as shown on the sticker sheet (Fig. 3).

### Reapplying a panel that has been removed

- 1 Rinse the panel in warm water; you may also use a mild detergent to clean it.
- 2 Let the panel air dry.
- 3 Reapply the panel to the white section on the front of the handle.

## Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 4).
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: There is a small gap between the brush head and the handle.

## Charging your Sonicare

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger (Fig. 5).

The battery charge indicator in the power on/off button flashes and the handle beeps twice to indicate that the toothbrush is charging.

### Battery low

If the battery charge of your Sonicare is low, you hear three beeps and the battery charge indicator flashes rapidly for 1 minute after the brushing cycle.

### Battery fully charged

When the Sonicare is fully charged, the battery charge indicator stops flashing and remains on for 10 seconds before it turns off.

To keep the battery fully charged at all times, you may keep your Sonicare on the charger when not in use. It takes at least 24 hours to fully charge the battery.

Note: When the battery charge is depleted, the handle will emit 2 sets of 5 beeps at the end of a brushing session or when a new brushing session is attempted. The new brushing session does not begin. The green LED will also flash fast until the handle goes to sleep.

## Bluetooth connectivity

Before you start the connectivity procedure, make sure your mobile device runs on Bluetooth 4.0 or higher.

Note: For the latest information on device compatibility, please visit the App Store or Google Play.

### Connecting your toothbrush

- 1 Download the Philips Sonicare For Kids app from the App Store or Google Play.
- 2 Confirm Bluetooth is switched on on your mobile device before you open the app. Check the settings on your mobile device to turn on Bluetooth.
- 3 To connect the toothbrush to the app, turn on the toothbrush by pressing either the power on/off button or the brushing mode button. Make sure the toothbrush is not on the charger when you turn it on.
- 4 Open the Philips Sonicare For Kids app.
- 5 Follow the app instructions displayed on your mobile device. The app automatically connects with the toothbrush. If the toothbrush does not connect, follow the troubleshooting steps provided in the app.
- 6 Keep your mobile device close to your toothbrush when you use your toothbrush. Make sure your mobile device is safe in a dry place.

## Using the appliance

### Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), where the teeth and the gumline meet. This is to ensure that the whole tooth surface is brushed (Fig. 6).
- 3 Press the power on/off button to turn on the Sonicare.
- 4 Apply light pressure to maximize Sonicare's effectiveness and let the Sonicare toothbrush do the brushing for you. Do not scrub (Fig. 7).
- 5 Maintain light pressure and gently brush the outer and inner surfaces of all top and bottom teeth while you keep the toothbrush at an angle such that the toothbrush bristles touch the area where the teeth and the gumline meet.

- 6 The inner and outer surfaces of the upper and lower back teeth are hard to reach. Spend time brushing these areas, apply light pressure and maintain the angle of the toothbrush bristles to the gumline (Fig. 8).
- 7 Brush the chewing surfaces of all the upper and lower teeth.
- 8 While you maintain the angle of the bristles to the gumline, gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between the teeth. Continue this motion throughout the brushing cycle (Fig. 9).

Note: To make sure your child brushes evenly throughout the mouth, divide the mouth into 4 sections using the KidPacer feature (see chapter "Features").

- 9 Begin brushing in section 1 (outside top teeth). A short series of tunes indicate it is time to move to section 2 (inside top teeth). At the next series of tunes, move to section 3 (outside bottom teeth). The next series of tunes indicate it is time to move to section 4 (inside bottom teeth). The toothbrush will automatically turn off at end of the brushing cycle (Fig. 10).

Your Sonicare is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

## Brushing modes

The Sonicare automatically starts in the default Gentle mode.

- 1 Press the brushing mode button to toggle between modes (Fig. 11). The green mode light indicates the selected mode.
- Gentle mode: Full power mode designed for older/independent brushers. When in Gentle mode, the brushing cycle is set at 2 minutes.

## 14 English

- Extra Gentle mode: Reduced power mode used for training and introducing younger children to the sonic experience. When in Extra Gentle mode, the brushing cycle is initially set at 1 minute and slowly progresses to 2 minutes over time to help young children reach the dentist recommended 2 minutes of brushing cycle.

To manually set the brushing cycle to 2 minutes in Extra Gentle mode, place the handle on the charger and hold down the mode button for 2 seconds. You hear a single beep to indicate the brushing cycle is set to 2 minutes.

To manually re-set the brushing cycle to 1 minute in Extra Gentle mode, place the handle on the charger and hold down the mode button for 2 seconds. You hear a double beep to indicate the brushing cycle is set to 1 minute.

## Features

### KidPacer

The Sonicare For Kids uses a short series of tunes to alert the brusher to move to the next brushing quadrant of the mouth and concludes the brushing cycle with a congratulatory tune before it automatically turns off. This ensures that the child's teeth get a thorough, overall clean (see "Brushing instructions").

Note: Instruct your child to brush until they hear the congratulatory tune and the toothbrush turns off.

### Easy-start

The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare. The Sonicare For Kids comes with the Easy-start feature activated.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

## Activating or deactivating the Easy-start feature

1 Place the handle on the charger.

To activate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated. To deactivate Easy-start: Press and hold the power on/off button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

## Cleaning

**Warning: Do not clean brush heads, handle, charger, charger cover or replaceable panel in the dishwasher.**

### Toothbrush handle

1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste (Fig. 12).

**Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.**

2 Wipe the entire surface of the handle with a damp cloth.

## Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 13).
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

## Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

## Storage

If you are not going to use the Sonicare for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

## Replacement

### Brush head

- Replace Sonicare brush heads every 3 months to achieve optimal results.
- Use only Sonicare For Kids replacement brush heads.

## Recycling

- This symbol means that this product shall not be disposed of with normal household waste (2012/19/EU) (Fig. 18).
- This symbol means that this product contains a built-in rechargeable battery which shall not be disposed of with normal household waste (Fig. 19) (2006/66/EC). We strongly advise you to take your product to an official collection point or a Philips service centre to have a professional remove the rechargeable battery.

- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

## Removing the rechargeable battery

**Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it.**

To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.

- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Insert a screwdriver into the slot located in the bottom of the handle and turn counterclockwise until you see a gap between the bottom cap and the handle (Fig. 14).
- 3 Insert the screwdriver into the gap and pry the bottom cap out of the handle (Fig. 15).
- 4 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 16).
- 5 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 17).

The rechargeable battery can now be recycled or disposed of and the rest of the product discarded appropriately.

## Guarantee and support

If you need information or support, please visit [www.philips.com/support](http://www.philips.com/support) or read the separate worldwide guarantee leaflet.

### Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads.
- Damage caused by use of unauthorised replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorised repair.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.
- Replaceable stickers
- Replaceable panels

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©2016 Koninklijke Philips N.V. (KPNV). All rights reserved.
Philips and the Philips shield are trademarks of KPNV.
Sonicare and the Sonicare logo are trademarks
of Philips Oral Healthcare, LLC and/or KPNV.