PHILIPS


FlexCare
Platinum
Connected

PHILIPS
sonicare

1

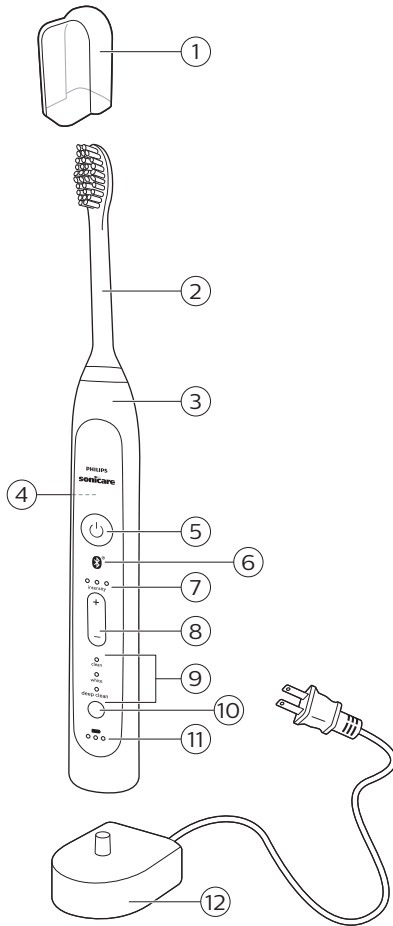— |

| —

— |

| —

English 6

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Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

Important safety information

Read this important information carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the wall socket.

Warning

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Always have the charger replaced with one of the original type in order to avoid a hazard.
- Do not use the charger outdoors or near heated surfaces.
- If the appliance is damaged in any way (brush head, toothbrush handle or charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see chapter 'Warranty and support').

Caution

- Do not clean the brush head, the handle and/or the charger in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use the Philips Sonicare.

- The Philips Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Philips Sonicare.
- This product has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the product and contact your doctor if you experience any discomfort or pain.
- The Philips Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Do not use other brush heads than the ones recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.

Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Your Philips Sonicare (Fig. 1)

- 1 Hygienic travel cap
- 2 Brush head with unique identifier
- 3 Handle
- 4 Sensors
- 5 Power on/off button
- 6 Bluetooth
- 7 Intensity level indicator
- 8 Intensity level button (+ and -)
- 9 Brushing modes

8 English

10 Mode button

11 Battery level indicator

12 Travel charger

- Not shown: Travel case

Note: The content of the box may vary based on the model purchased.

Preparing for use

Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle (Fig. 2).
- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: Philips Sonicare brush heads are imprinted with icons to easily identify your brush head (Fig. 4).

Charging your Philips Sonicare

- 1 Put the plug of the charger in a wall socket.
- 2 Place the handle on the charger.

The flashing light of the battery level indicator indicates that the toothbrush is charging (Fig. 5).

Note: It takes at least 24 hours to fully charge the battery.

Using the appliance

Bluetooth connectivity

Before you start the connectivity procedure, make sure your mobile device runs on Bluetooth 4.0 or higher.

Note: For the latest information on device compatibility, please visit the App Store or Google Play.

Connecting your toothbrush

- 1 When used in connected mode, the Philips Sonicare app gives you the opportunity to receive real-time feedback on brushing, store and share brushing data with your dentist (if you so choose), and receive personalized coaching. If you wish to use the toothbrush in unconnected mode, the toothbrush will still function. To use the toothbrush in connected mode, please follow the instructions below.
- 2 Download the Philips Sonicare app from the Apple App Store or Google Play Store.
- 3 Confirm Bluetooth is switched on on your mobile device before you open the app. Check the settings on your mobile device to switch on Bluetooth.
- 4 To connect the toothbrush to the app, switch on the toothbrush by pressing either the power on/off button or the brushing mode button. Make sure the toothbrush is not on the charger when you switch it on.
- 5 Open the Sonicare app.
- 6 Follow the Sonicare app instructions displayed on your mobile device. The app automatically connects with the toothbrush. If the toothbrush does not connect, follow the troubleshooting steps provided in the app.

Note: For optimal data capture, keep your mobile device close during brushing.

Brushing instructions

- 1 Wet the bristles and apply a small amount of toothpaste.
- 2 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline (Fig. 6).

Note: Keep the center of the brush in contact with the teeth at all times.

- 3 Press the power on/off button to turn on the Philips Sonicare.

10 English

- 4 Apply light pressure to maximize Philips Sonicare's effectiveness and let the Philips Sonicare toothbrush do the brushing for you. Do not scrub.
- 5 Gently move the brush head slowly across the teeth in a small back and forth motion so the longer bristles reach between your teeth. Continue this motion throughout the brushing cycle.
- 6 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth (Fig. 7).

Note: To make sure you brush evenly throughout the mouth, divide your mouth into 6 sections using the BrushPacer feature (see chapter 'Features').

- 7 Follow the Sonicare app's brushing guidance combined with the BrushPacer and begin brushing in section 1 (upper right teeth) and brush for 20 seconds before you move to section 2 (upper front teeth). Continue doing this for section 3 till 6 to ensure even brushing in all 6 segments of the mouth (Fig. 8).

Note: When using the White mode, two additional 20 second intervals can be used for polishing your front teeth.

- 8 The Philips Sonicare app helps to address any missed spots by guiding you to spend 20 seconds in the areas you missed during brushing.

Your Philips Sonicare toothbrush is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Note: Ensure that braces or dental restorations are properly adhered and not compromised. If not, braces or dental restorations may be damaged when they come into contact with the brush head.

Note: When the Philips Sonicare toothbrush is used in clinical studies, the handle should be fully charged and the Easy-start feature deactivated.

Brushing modes

Note: When the toothbrush is switched on, you cannot toggle between modes.

- 1 Prior to turning on the Sonicare, press the mode button to toggle between modes.
 - The selected brushing mode illuminates.
 - Once a brushing mode has been selected, it remains active until another brushing mode is selected.

Clean mode

Standard mode for superior teeth cleaning.

White mode

2 minutes of alternating between White and Massage, with an additional 40 seconds to polish your visible front teeth.

Deep Clean

Provides an invigorating of deep massage to deliver an exceptionally clean experience.

Intensity levels per mode

- Low: 1 LED indicator light.
- Medium: 2 LED indicator lights.
- High: 3 LED indicator lights.

Note: to adjust the intensity level to the next higher or lower level, you can press the + button to increase intensity or the - button to decrease intensity at any time during the brushing cycle.

Pressure Sensor

A change in vibration (and a slight change in sound) alerts you when you apply too much pressure while brushing. See chapter 'Features'.

Battery status (when handle is not on charger)

- 1 yellow LED indicates a low charge.
- 3 green LEDs indicate a full charge.

Note: If you hear two sets of five beeps and the yellow LED of the battery level indicator flashes, your Sonicare will turn off until placed on the charger, where it will start to recharge.

Note: Your Sonicare will not be able to run a brushing routine until it is fully recharged.

Note: To keep the battery fully charged at all times, you may keep your Sonicare on the charger when not in use.

Features

Pressure Sensor

This Philips Sonicare model comes with the Pressure Sensor feature activated. The Pressure Sensor works in all 3 modes and at all 3 intensity levels.

Activating or deactivating the Pressure Sensor feature

- 1 Attach the brush head to the handle.
- 2 Place the handle in the plugged-in charger.
 - To deactivate Pressure Sensor:
Press and hold the power on/off and mode button for 2 seconds. You hear 1 beep to indicate that the Pressure Sensor feature has been deactivated.
 - To activate Pressure Sensor:
Press and hold the power on/off and mode button for 2 seconds. You hear 2 beeps to indicate that the Pressure Sensor feature has been activated.

BrushPacer

- The BrushPacer is an interval timer that has a short beep and pause to remind you to brush the 6 sections of your mouth. Depending on the brushing mode or brushing routine you have selected, the BrushPacer beeps at different intervals during the brushing cycle. See the 'Brushing instructions' sections (Fig. 8).
- To deactivate or reactivate the BrushPacer feature:
- 1 Place the handle with a brush head attached in the plugged-in charger.
 - 2 Press and hold the mode button until:
 - a You hear 1 beep to indicate the BrushPacer has been deactivated.
 - b You hear 2 beeps to indicate the BrushPacer has been activated.

Smartimer

The Smartimer indicates that your brushing cycle is complete by automatically turning off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Easy-start

This Philips Sonicare comes with the Easy-start deactivated.

The Easy-start feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly advance through the Easy-start ramp-up cycle.

Activating or deactivating the Easy-start feature

- 1 Attach the brush head to the handle.
 - 2 Place the handle in the plugged-in charger.
- To activate Easy-start:

Press and hold the power on/off and mode button for 2 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

- To deactivate Easy-start:

Press and hold the power on/off and mode button for 2 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.

Note: Using the Easy-start feature beyond the initial ramp-up period is not recommended and reduces Philips Sonicare effectiveness in removing plaque.

Note: Each of the first 14 brushings must be at least 1 minute in length to properly move through the Easy-start ramp-up cycle.

Cleaning

Warning: Do not clean the brush head, handle or the travel charger in the dishwasher.

Toothbrush handle

- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste (Fig. 9).

Caution: Do not push on the rubber seal on the metal shaft with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.

Brush head

- 1 Rinse the brush head and bristles after each use (Fig. 10).
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week. Rinse the travel cap as often as needed.

Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the product for a long time, unplug it, clean it (see chapter 'Cleaning') and store it in a cool and dry place away from direct sunlight.

Replacement

Brush head

- Replace Philips Sonicare brush heads every 3 months to achieve optimal results.

Note: The Sonicare App will tell you when it is time to replace your brush head. (See in-app Help section)

- Use only Philips Sonicare replacement brush heads.

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

- This product contains a built-in rechargeable battery which shall not be disposed of with normal household waste. Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

Removing the rechargeable battery

Warning: PLEASE NOTE THAT THIS PROCESS IS NOT REVERSIBLE.

- To remove the rechargeable battery, you need a flat-head (standard) screwdriver. Observe basic safety precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.
- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer switch on the Philips Sonicare.
 - 2 Insert a screwdriver into the slot located at the bottom of the handle and turn counterclockwise to release the bottom cap (Fig. 11).
 - 3 Hold the handle upside down and push down on the shaft to release the internal components of the handle (Fig. 12).
 - 4 Insert the screwdriver under the circuit board, next to the battery connections, and twist to break the connections. Remove the circuit board and pry the battery from the plastic carrier (Fig. 13).

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Warranty and support

If you need information or support, please visit

www.philips.com/support or read the international warranty leaflet.

Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Hereby, Philips Consumer Lifestyle BV declares that the radio equipment type Bluetooth LE is in compliance with Directive 2014/53/EU.

The full text of the EU declaration of conformity is available at the following internet address: **www.philips.com/support**

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©2018 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.