



PHILIPS

Colgate

SonicPro 50

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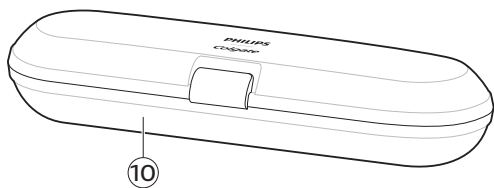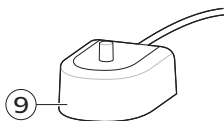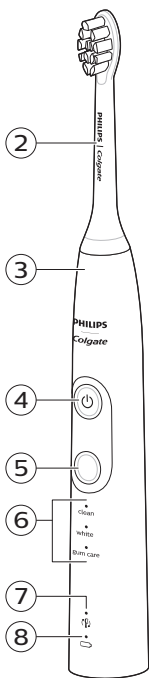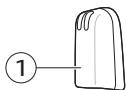Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger away from water. Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not immerse the charger in water or any other liquid.
- After cleaning, make sure the charger is completely dry before you connect it to the wall socket.

Warning

- The mains cord cannot be replaced. If the mains cord is damaged, discard the charger.
- Always have the charger replaced with one of the original type in order to avoid a hazard.
- Do not use the charger outdoors or near heated surfaces.
- If the appliance is damaged in any way (brush head, toothbrush handle or charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Consumer Care Center in your country (see chapter 'Warranty and support').
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.

Caution

- Do not clean the brush head, travel cap, handle, charger or travel case in the dishwasher.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use. Also consult your dentist if you experience discomfort or pain when you use the Philips Colgate toothbrush.
- The Philips Colgate toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device prior to use.
- If you have medical concerns, consult your doctor before you use the Philips Colgate toothbrush.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear. Do not use other brush heads than the ones recommended by the manufacturer.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristle tufts dislodging.

- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of this product's plastic components.

Battery safety instructions

- Only use this product for its intended purpose and follow the general and battery safety instructions as described in this user manual. Any misuse can cause electric shock, burns, fire and other hazards or injuries.
- Only use the charger type HX6100/HX6160 provided with the appliance to charge the battery.
- Do not use an extension cord with this product.
- Charge, use and store the product at a temperature between 0 °C and 40 °C.
- Always unplug the product when it is fully charged.
- Do not burn products and their batteries and do not expose them to direct sunlight or to high temperatures (e.g. in hot cars or near hot stoves). Batteries may explode if overheated.
- If the product becomes abnormally hot, gives off an abnormal smell, changes color or if charging takes much longer than usual, stop using and charging the product and contact Philips directly via a service center or warranty line at www.philips.com/support, or toll free at 1-800-243-3050.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- This product contains a rechargeable battery that is non-replaceable. Do not open the product to replace the rechargeable battery during the life of the device.

- When you handle batteries, make sure that your hands, the product and the batteries are dry.
- To prevent batteries from heating up or releasing toxic or hazardous substances, do not modify, pierce or damage products and batteries and do not disassemble, short-circuit, overcharge or reverse charge batteries during the life of the device.
- To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.
- If disposable batteries are damaged or leaking, avoid contact with the skin or eye. If this occurs, immediately rinse well with water and seek medical care.

Electromagnetic fields (EMF)

This Philips Colgate toothbrush complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Introduction

Congratulations on your purchase and welcome to Philips Colgate! To fully benefit from the support that Philips Colgate offers, register your product at **www.philips.com/welcome**.

Your Philips Colgate toothbrush (Fig. 1)

- 1 Hygienic travel cap
- 2 Smart brush head(s)*
- 3 Handle
- 4 Power on/off button
- 5 Mode/intensity button
- 6 Modes
- 7 Brush head replacement reminder light
- 8 Battery indication
- 9 Charging base
- 10 Travel case (specific types only)*

*Note: The content of the box may vary based on the model purchased.

Brush heads

Your Philips Colgate toothbrush comes with one or more brush heads which are designed to deliver superior results for your oral care needs.

BrushSync Technology

BrushSync technology enables your brush heads to communicate with your handle using a microchip. This enables the brush head replacement reminder feature.

Brushing modes

To toggle between modes press the mode button located below the on/off button.

Your power toothbrush is equipped with 3 brushing modes, depending on your model. The following modes are available:

Clean mode

Clean is a 2-minute mode. It gives you a thorough and complete clean.

White mode

White is a 2-minute 30-seconds mode.

Gum Care mode

Gum Care is a 3-minute mode.

Battery status and charging

Charging

- 1 Plug the charger into an electrical outlet.
- 2 Place the handle on the charger.
 - The handle will emit two beeps when charging begins.
 - The battery indicator light will flash when the toothbrush is charging.
 - When the handle is fully charged on the charger, the battery indicator light will turn solid green for 30 seconds and switch off.

Note: Your toothbrush comes pre-charged for first use. After first use, charge for at least 24 hours.


Battery status (when handle is not on charger)

When removing the handle from the charger, the battery indicator light at the bottom of the handle will indicate the status of the battery.

- Solid green LED: full battery
- Flashing green LED: medium battery
- Flashing amber LED and three beeps: low battery
- Flashing amber LED and two sets of five beeps: No brushing sessions left (charge toothbrush)


Using your Philips Colgate toothbrush

Brushing instructions

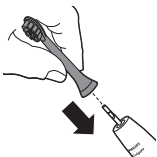- 1 Firmly press the brush head down onto the metal shaft of the handle until it stops. The brush head bristles should be facing the same direction as the front of the handle.


Note: It is normal to see a slight gap between the brush head and the handle.

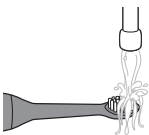- 2 Wet the bristles.


- 3 Apply a small amount of toothpaste on the bristles.

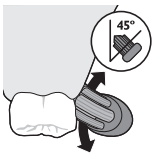- 4 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing firmly to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

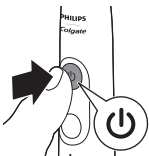- 5 Press the power on/off button to turn on the Philips Colgate toothbrush.


- 6** Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth.

Note: The bristles should slightly flare. Do not scrub.


- 7** To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

Note: A change in vibration (and a slight change in sound) alerts you when you apply too much pressure while brushing.

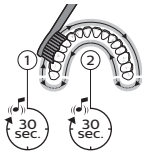- 8** To make sure you brush evenly throughout the mouth, divide the mouth into 4 sections using the Quadpacer feature. The Quadpacer feature beeps every 30 seconds over 2 minutes.

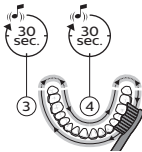- 9** Brush for the full 2 minutes following the Quadpacer.

In Clean mode, after 2 minutes the handle will automatically shut off.

In White mode, an additional 30 seconds is provided for you to brighten and polish your front teeth.

In Gum care mode, after 2 minutes the handle will continue to brush for an additional minute so you can gently stimulate and massage your gums.

- After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs.

Note: Your Philips Colgate toothbrush is safe to use on:

- Braces (brush heads wear out sooner when used on braces)
- Dental restorations (fillings, crowns, veneers)

Note: When the Philips Colgate toothbrush is used in clinical studies, it must be used in the Clean mode (for plaque and gum health claims) or in White mode (for whitening claims) at high intensity with EasyStart turned off.

Features

- EasyStart
- Brush head replacement reminder
- Pressure sensor
- Quadpacer


EasyStart

The EasyStart feature gently increases the power of the handle over the first 14 brushings to help you get used to the brushing power with the Philips Colgate toothbrush.

The EasyStart feature works in all modes.

Note: The EasyStart feature comes activated with your product.

To deactivate EasyStart see 'Activating or deactivating features'.

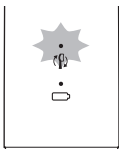Brush head replacement reminder

Your Philips Colgate toothbrush is equipped with BrushSync technology that tracks the wear of your brush head.

- When attaching a new smart brush head for the first time the brush head replacement reminder light blinks green three times. This confirms you have a Philips brush head with BrushSync technology.
- Over time, based on the overall pressure you apply while brushing and the total time you have brushed with your brush head, the handle will track brush head wear in order to determine the total time you have brushed with your brush head. This feature gives you the guarantee for the best cleaning and care of your teeth.
- When the brush head replacement reminder light lights up amber, you should replace your brush head.

Note: The brush head replacement reminder comes activated with your product. To deactivate the brush head replacement reminder see 'Activating or deactivating features'.

Pressure sensor

Your Philips Colgate toothbrush measures the pressure you apply while brushing to protect your gums and teeth from damage. If you apply too much pressure, the handle will change in vibration until you reduce the pressure.

Note: The pressure sensor comes activated with your product. To deactivate the pressure sensor see 'Activating or deactivating features'.

Quadpacer

The Quadpacer is an interval timer that has a short beep and pause to remind you to brush the different sections of your mouth. Depending on the brushing mode you have selected (see 'Brushing modes'), the Quadpacer beeps at different intervals during the brushing cycle.

- In Clean mode, the Quadpacer beeps every 30 seconds over the 2 minute brushing cycle.
- In White mode, the Quadpacer beeps every 30 seconds over the first 2 minutes and beeps every 15 seconds for the final 30 seconds of the brushing cycle.
- In Gum care mode, the Quadpacer beeps every 30 seconds over first 2 minutes and beeps every 15 seconds for the final minute of the brushing cycle.

Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder
- Pressure sensor

To activate or deactivate these features, follow the instructions below:

Step 1: Place the handle on the charger.

Step 2: Press and hold the power button for:

Easy start

Brush head
replacement
reminder

Pressure sensor

Up to 3
seconds

Up to 5
seconds

Up to 7 seconds

Step 3: Let go power button when you hear

1 beep

2 beeps

3 beeps

If you see the brush head replacement light and battery indication light blink green 2 times, and hear 3 tones from low to high, then the feature has been activated.

OR

If you see the brush head replacement light and battery indication light blink amber once, and hear 3 tones from high to low, then the feature has been deactivated.

Cleaning

Warning: Do not clean the brush head, travel cap, handle, charger or travel case in the dishwasher.

Toothbrush handle


- 1 Remove the brush head from the handle and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste.

Caution: Do not push on the rubber seal at the top of the handle with sharp objects, as this may cause damage.

- 2 Wipe the entire surface of the handle with a damp cloth.


Note: Do not use isopropyl rubbing alcohol, vinegar or bleach to clean handle as this may cause discoloration.

Brush head

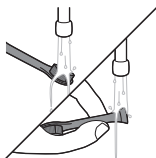- 1 Remove the brush head from the handle and rinse it thoroughly.
- 2 Rinse the brush head connection with warm water at least once a week.

Charger

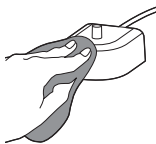- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth. Make sure you remove any toothpaste or other residue from the charger surfaces.

Storage

If you are not going to use the product for an extended period of time, unplug the charger from the electrical outlet, clean it and store it in a cool and dry place away from direct sunlight.

Locating the model number

Look on the bottom of the Philips Colgate toothbrush handle for the model number (HX684x).

Recycling

- Do not throw away the product with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- This product contains a built-in rechargeable battery which shall not be disposed of with normal household waste. Please take your product to an official collection point or a Philips service center to have a professional remove the rechargeable battery.
- Follow your country's rules for the separate collection of electrical and electronic products and rechargeable batteries. Correct disposal helps prevent negative consequences for the environment and human health.

Removing the rechargeable battery


Warning: Only remove the rechargeable battery when you discard the appliance. Make sure the battery is completely empty when you remove it. To remove the rechargeable battery, you need a towel or cloth, a hammer and a flat-head (standard) screwdriver. Observe basic safety

precautions when you follow the procedure outlined below. Be sure to protect your eyes, hands, fingers, and the surface on which you work.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Colgate toothbrush and let it run until it stops. Repeat this step until you can no longer turn on the Philips Colgate toothbrush.
- 2 Remove and discard the brush head. Cover the entire handle with a towel or cloth.


- 3 Hold the top of the handle with one hand and strike the handle housing 0.5 inch above the bottom end. Strike firmly with a hammer on all 4 sides to eject the end cap.

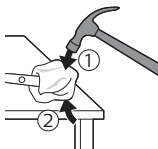
Note: You may have to hit on the end several times to break the internal snap connections.

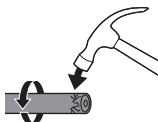
- 4 Remove the end cap from the toothbrush handle. If the end cap does not release easily from the housing, repeat step 3 until the end cap is released.

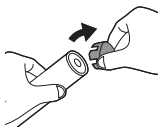
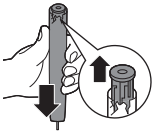
- 5 Holding the handle upside down, press the shaft down on a hard surface. If the internal components do not easily release from the housing, repeat step 3 until the internal components are released.

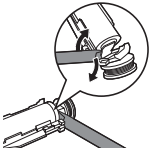- 6 Wedge the screwdriver between the battery and the black frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the black frame.

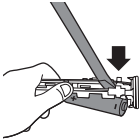- 7 Insert the screwdriver between the bottom of the battery and the black frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

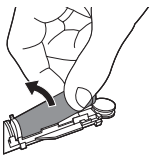
- 8 Grab the battery and pull it away from the internal components to break the second metal battery tab.

Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.

- 9 Cover the battery contacts with tape to prevent any electrical short from residual battery charge.

The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

Note: Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.

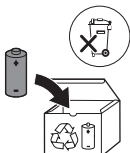Warranty and support

If you need information or support, please visit www.philips.com/support or read the international warranty leaflet.

WARRANTY EXCLUSIONS

What is not covered under warranty:

- Brush heads.
- Damage caused by use of unauthorized replacement parts or unauthorized brush heads.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

Colgate and the Colgate Logo are the registered trademarks of Colgate-Palmolive Company.

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©2019 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV.