

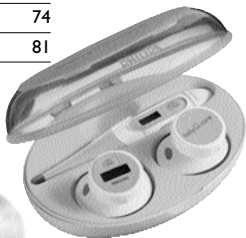
# SCH530

<b>EN</b>	User manual	4
<b>FR</b>	Mode d'emploi	11
<b>ES</b>	Instrucciones de manejo	18
<b>DE</b>	Bedienungsanleitung	25
<b>NL</b>	Gebruiksaanwijzing	32
<b>IT</b>	Istruzione per l'uso	39
<b>PR</b>	Manual de utilização	46
<b>GR</b>	Οδηγίες χρήσεως	53
<b>SV</b>	Bruksanvisning	60
<b>DK</b>	Vejledning	67
<b>NO</b>	Bruksanvisning	74
<b>SU</b>	Käyttöohje	81
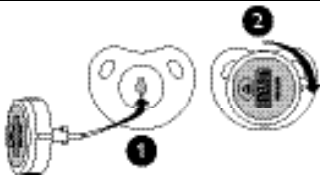
# PHILIPS
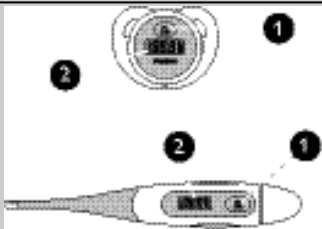
### **Soother thermometer**

The most convenient way to measure your baby's temperature is to use the soother thermometer. This is also the most comfortable method for baby. Don't use for 30 minutes after the baby has had hot or cold drinks as this causes an inaccurate measurement.

- Push the on/off button (A-1) to switch on the thermometer. (The display (A-2) shows an 'L' and a flashing '°C'.)
- Place the soother in the baby's mouth and leave it there until the thermometer beeps four times.
- The display now shows the measured temperature and the °C sign lights continuously.
- Remove the soother from the baby's mouth.
- Push the on/off button to switch off the thermometer. (Doing this extends battery power. The thermometer switches off automatically after about 10 minutes.)

### **Checking the last measured temperature**

If you wish to compare temperatures then you can check to see what the last measured temperature was. This can only be done once before a new measurement. Push and hold the on/off button for about 3 seconds. The display shows the last measured temperature and an 'M' for memory. When you release the button this temperature is erased.

### **Battery low indication**

This soother thermometer has a long-life battery which should last for about 2 to 3 years. If the display shows a , the battery is low and the thermometer has reached the end of its useful life. Do not use it any more.

Before disposing of the thermometer remove the battery and dispose of it in an environmentally friendly manner. (To gain access to the battery unclip the small cover.)

### **Cleaning the soother**

Nipple, shield and cover can be sterilised\* or boiled. Never sterilise nor boil the measurement unit. Before cleaning always remove the measurement unit from the nipple and shield as follows (see figure B):

- Turn the measurement unit  $90^\circ$  anticlockwise with respect to the shield.
- Pull the measurement unit straight out of the nipple and shield. (Be careful not to bend the tip of the measurement probe.)

After boiling the nipple, shield and cover, drain and push out any remaining water and allow to dry.

To clean the measurement unit wipe it with a damp cloth.

Reassemble the soother as follows (see figure C):

- Align the blocks of the measurement unit with the matching slots in the back of the shield.
- Insert the measurement unit into the shield ensuring that the tip of the probe goes into the pocket inside the nipple.
- Turn the measurement unit  $90^\circ$  clockwise to lock into position.

\* Sterilising means destroying harmful bacteria at a temperature of  $95^\circ\text{C}$  for more than 3 minutes. It does NOT mean sterilisation according to clinical standards. Philips Baby Bottle Steriliser is suitable for sterilising the nipple, shield and cover.

**Digital flexible-tip thermometer**

The soother thermometer is as accurate as the digital flexible-tip thermometer and is more convenient to use. However, because rectal measurement is inherently more accurate than oral measurement, the digital flexible-tip thermometer gives you an opportunity to reassure yourself.

Remove the digital flexible-tip thermometer from its storage cover. If available, use special plastic hygiene covers to cover the tip of the thermometer (these are available from most pharmacies).

For a rectal measurement of a baby's temperature proceed as follows:

- Lubricate the tip of the thermometer (or hygiene cover) with a water-soluble jelly. Do not use a petroleum jelly.
- Push the on/off button (A-1) to switch on the thermometer. (The display (A-2) shows an 'L' and a flashing '°C'.)
- Lie the baby on its stomach or its side.
- With one hand separate the baby's buttocks. With the other hand gently slide the tip of the thermometer no more than 12mm (1/2 inch) into the rectum. If you detect any resistance stop immediately.
- Hold the baby firmly so that it cannot turn while the thermometer is inserted.
- Keep the thermometer in place until it beeps four times.
- Remove the thermometer and read the display.
- Push the on/off button to switch off the thermometer. (Doing this extends battery power. The thermometer switches off automatically after about 10 minutes).

**Checking the last measured temperature**

If you wish to compare temperatures then you can check to see what the last measured temperature was. This can only be done once before a new measurement. Push and hold the on/off button for about 3 seconds. The display shows the last measured temperature and an 'M' for memory. When you release the button this temperature is erased.

### Measuring the temperature of older children

For older children and adults the digital flexible-tip thermometer can also be used to measure temperature in the mouth (oral) or under the arm (the least reliable method). If the thermometer has been used rectally without hygienic covers, we do not recommend using it for oral measurements.

### Cleaning

Clean the thermometer after each use. Never sterilise nor boil the thermometer.

For home use, wipe the thermometer with a soft damp cloth or wash in water. The thermometer is water resistant.

Do not wash in the dishwasher.

### Replacing the battery

The thermometer has a long-life battery which should last for about 2 to 3 years. If the display shows a , the battery is low. Replace the battery as follows before using the thermometer:

- Remove the battery cover.
- Use an insulated probe to remove the battery (a toothpick for example).
- Insert the new battery (LR41) with the positive side (+) facing upwards.
- Secure the cover back in place.

### Warning

Swallowing a button battery can be extremely dangerous.

Dispose of battery in an environmentally friendly way.

Do not dispose of in fire or it may explode.

Never attempt to recharge the battery.

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Trademarks are the property of Koninklijke Philips Electronics N.V.
2005 © Koninklijke Philips Electronics N.V. All rights reserved.
**Philips DAP B.V.**